AT51047
INSTRUKCJA UŻYTKOWANIA CHODZIK
USER MANUAL ROLLATOR
BEDIENUNGSANLEITUNG ROLLATOR
NÁVOD K POUŽITÍ CHODÍTKA
NÁVOD NA POUŽITIE CHODÍTKO
MANUEL D'INSTRUCTION DÉAMBULATEUR
HANDLEIDING ROLLATOR
MANUAL DE INSTRUCCIONES ANDADOR
MANUALE DI ISTRUZIONI DEL DEAMBULATORE
BRUKSANVISNING RULLATOR

PI

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EN

Thank you for purchasing our product. Please read the user manual carefully before use. Do not operate the product without first understanding and becoming familiar with this manual. If you do not understand the warnings, precautions, or recommendations, consult a healthcare professional or the seller to avoid product damage or personal injury. **WARNING:** Check all components for damage caused during transport. If any damage is detected, DO NOT use the product. For more information, contact the seller.

INDICATIONS

This product is designed to assist individuals with limited walking ability. The rehabilitation walker is intended for single-person use. It is suitable for the following purposes: The rehabilitation walker allows patients to practice balance, which is crucial in the rehabilitation process after injuries or neurological disorders. Patients can use the rehabilitation walker to train correct posture, which is especially important for those with postural defects or stability issues. Exercises performed with the rehabilitation walker can help improve motor coordination and movement control, essential for individuals with neurological disorders.

The rehabilitation walker can aid in the rehabilitation process after lower limb injuries, helping patients gradually regain mobility.

For some patients, the rehabilitation walker can increase independence in daily life, enabling autonomous movement.

CONTRAINDICATIONS

Physical or mental limitations (e.g., impaired vision) that prevent the safe handling of the product.

Individuals with severe spinal stability issues or spinal injuries may not be suitable candidates for the rehabilitation walker, as it may increase the risk of further injuries. Using the rehabilitation walker may be unsafe for individuals with uncontrolled seizures or epilepsy, posing risks to both the user and others. Individuals with very limited muscle strength may find using the rehabilitation walker challenging or ineffective. Alternative rehabilitation approaches, such as manual therapy or strength-enhancing exercises, may be necessary.

For patients with severe balance disorders, using the rehabilitation walker may increase the risk of falls and further injuries.

FEATURES

Easy height adjustment

Shopping bag

Lightweight frame

Comfortable seat

Parking brakes

Made of carbon fiber

TECHNICAL SPECIFICATION

Overall length: 65 cm

Overall width (unfolded): 65 cm

Overall width (folded): 25 cm

Height when folded (lowest): 84 cm

Folded dimensions (L × W × H): 65 × 25 × 84 cm

Adjustable handle height (min–max): 78–98 cm (7 levels)

Distance between handles: outer: 54 cm inner: 43 cm

Seat width: 46 cm

Seat depth: 20.5 cm

Seat height (from ground): 63.5 cm

Front wheels: 200 × 36 mm (TPR)

Rear wheels: 200 × 36 mm (TPR)

Total weight (without bag and backrest): 5.57 kg

Backrest weight: 0.165 kg

Mesh bag weight: 0.28 kg

Maximum user weight: 150 kg

Maximum bag load: 5 kg

THIS SYMBOL INDICATES THE MAXIMUM USER WEIGHT

COMPONENTS

Rehabilitation Walker, Backrest strap, User manual

IMPORTANT:

In the event of a "serious incident" related to the product that directly or indirectly led to, could have led to, or may lead to any of the following: a) Death of the patient, user, or another person, or

b) Temporary or permanent deterioration of health of the patient, user, or another person, or

c) Serious threat to public health,

The "serious incident" must be reported to the manufacturer and the competent authority in the user's or patient's country of residence.

WARNING: If pain, allergic reactions, or other concerning and unclear symptoms occur while using this medical device, consult a healthcare professional.

WARNING:

Ensure no stones or other objects are blocking the walker's wheels.

Before using the walker as support for standing up, make sure both brake levers are in the parking position. Once the user is standing comfortably and stably, the parking brakes can be released.

Keep both hands on the handles while walking.

Remember that some parts of the walker may become very hot or cold due to ambient temperature, sunlight, heating devices, etc.

Exercise caution when touching them.

Before each use, check that the handles are securely attached.

Do not hold the handles with wet hands. In such cases, you may lose grip and balance.

Do not drive the walker onto the road. Always stay on the sidewalk.

Be careful not to trap clothing or other items in the frame of the walker when folding/unfolding it.

WARNING: It is prohibited to use the product in any way other than its intended purpose.

WARNING: The manufacturer is not responsible for damages caused by neglecting cleaning and maintenance rules or failing to follow the recommendations in this user manual.

TARGET PATIENT GROUP

Individuals with diseases, dysfunctions, or injuries for which this product is designed for treatment/rehabilitation or compensation (refer to the section on product purpose in this user manual). The product can be purchased independently by the user or based on the recommendation of a doctor, therapist, or other specialist. Regardless of whether the product is purchased independently or upon recommendation, consideration should be given to the available sizes/necessary features, product variants, indications and contraindications for use, as well as information provided by the manufacturer.

OPERATING INSTRUCTIONS

Unfolding:

Remove the rubber cover from the handles.

Pull the handles outward.


Push down the outer seat rails until you hear a distinct locking click.


Backrest Assembly:

Follow the instructions included in the manual for proper installation.

①

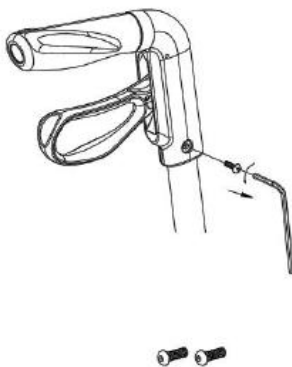②

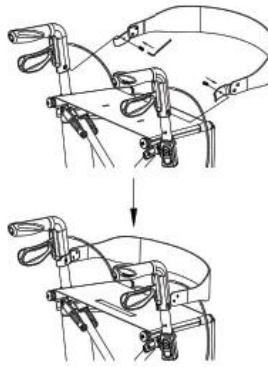Folding:

Lift the seat using the nylon strap in the middle to unlock the mechanism.

Push the handles together to ensure complete folding.

Height Adjustment:

Press the lock button.

Adjust the handles to the desired height.

WARNING: Both handles should be set to the same height.

**Braking System:**

The walker is equipped with brakes operated by two brake levers, offering two modes: stop brakes and parking brakes.

Stop Brakes: Pull both brake levers upward. Release to stop braking.

Parking Brakes: Push both brake levers down to the locked position. Always engage both parking brakes simultaneously. To disengage, pull both brake levers upward.

Sitting, Standing, and Walking:

WARNING: Ensure both brake levers are in the parking position before sitting.

Engage the parking brakes by pressing both levers down.

Use the handles for support, turn around, and sit on the seat.

To stand, push on the handles for support, turn, and disengage the parking brakes.

Release the brake levers and step forward, keeping the walker close.

MAINTENANCE AND CLEANING

Clean all surfaces with a cloth and clean water; if necessary, use a small amount of mild detergent. Dry with a soft cloth.

Avoid using harsh detergents or abrasive materials to prevent scratches.

Before each use, ensure the walker is in good condition.

Regularly inspect:

Wheels (freedom of rotation and condition)

Frame (stability and connections)

Handles (secure height adjustment and support)

Seat (wear or damage)

STORAGE

Store the product in a cool, dry place, away from moisture and direct sunlight.

DISPOSAL

Dispose of the product as regular household waste following local recycling guidelines. Protect the environment.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp. z o.o., 03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR MEDICAL Sp. z o.o.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

Date of issue of the manual: 18.11.2024
v1-18.11.2024 v2-02.01.2026 v3-23.03.2026

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