AT51031

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USER MANUAL ROLLATOR

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The outdoor rollator is equipped with a parking brake and an upholstered seat with backrest. The height of the handle is adjustable and the rollator can be folded for transport or storage.

Before assembling or using the rollator for the first time, read the operating instructions carefully and keep them in a safe place.

DESTINY

The device serves as a walking aid, support and carrying smaller items both inside and outside. It is intended to be used by a user with reduced mobility.

In order to guarantee the safety of the user of the baby rollator, the user should have adequate physical, cognitive and visual abilities. If in doubt, please contact the retailer / specialist who recommended the use of the rollator.

It is very important that the user, before using the device for the first time, carefully reads the safety instructions contained in this manual. The instructions listed here are intended to prevent possible injury to the user or a third person as well as damage to the rollator or other objects within it.

TIP:

It is impossible to stand on individual parts of the frame. The frame of the device is not designed for this purpose.

Contraindications

physical or mental limitations (e.g. visual impairment) that prevent the safe handling of the product.

SECURITY

The rollator can be used by a person weighing a maximum of 136 kg.

- (1) Please make sure that the brake of an unattended rollator is always locked / on.
- (2) Before use, please make sure that all wheels, including moving wheels, are securely fastened.
- (3) Please make sure the brakes are functioning properly.
- (4) Before using for the first time, please check the height of the handle and its correct fit.
- (5) Please make sure the rollator is fully unfolded.
- (6) When unfolding the rollator, the frame should be set in the correct position.
- (7) To avoid pinching your fingers, please stay away from the folding mechanism.
- (8) In order to adapt the rollator to your needs, please seek advice from qualified personnel.
- (9) Exercise care when lifting or moving the rollator.

CAUTION:

When a product-related „serious incident“ occurs, which directly or indirectly led, or could lead to, any of the following events:

- a) death of a patient, user or other person or
- b) temporary or permanent deterioration in the health of a patient, user or other person, or
- c) a serious threat to public health

the above „serious incident“ must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient is resident. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

CAUTION:

In the event of pain, allergic reactions or other disturbing, unclear for the user symptoms related to the use of a medical device, consult a healthcare professional.

FOLDING AND STORING THE PADDLE

1. Pull the strap on the center of the saddle upwards to fold down the rollator.

IT IS FORBIDDEN:

- (1) Using the rollator as a means of transporting people - the appliance is not designed as a wheelchair.
- (2) Do not sit on the rollator without the brake locked - this could set the rollator in motion and could injure the user or a third party.
- (3) Do not sit on a rollator on a slope without activating the parking brake.
- (4) Do not sit on the saddle with a flourish - please sit carefully and in a controlled manner, ask for help if necessary.
- (5) Do not hang any objects on the brake rod or on the handles, this may impair the operation of the brakes and / or the stability of the rollator.
- (6) Please do not step on the rollator or use it as a support for the feet.

INSTALLATION INSTRUCTIONS

Please make sure that the rollator has not been damaged during transport. In case of visible damage, please do not use the rollator. In this case, please contact your dealer to have the rollator repaired or returned.

1. Remove the packaging. The rollator is folded up.
2. Unfold the rollator, press the saddle down on both sides.

NOTE: Both sides must be locked into position with an audible „click“.

3. Hang the bag.

HEIGHT ADJUSTMENT

Adequate customization

Before using the rollator for the first time, please make sure that it is fully extended and that the frame is properly locked and stiff. It is important to adjust the height of the rollator correctly in a comfortable position.

1. To do this, lower your hands and set the handles of the rollator at the level of your wrists.
2. When holding the handles, your elbows should be bent at an angle of 15 ° -20 °.
3. If you have to bend down or raise your arms to reach the rollator, please correct the height.
4. In doing so, remove the thumbscrew and adjust the handle to the correct height. NOTE: The adjusting screw should be screwed in ONLY after you hear a click. Screwing the screw beyond the holes not only scratches the frame, but also leads to mechanical damage, which is not covered by the manufacturer's warranty. Please repeat this on the other side of the rollator and check before use that both handles are locked at the same height and correctly.

TIP: To determine the correct height of the rollator, please consult your physician or specialist dealer.

BRAKES SETTING AND MAINTENANCE

Service brake (when moving): Pull the brake lever upwards (fig. 2).

Parking brake (when parking): Press the brake lever downwards (fig. 1).

Brake adjustment: To increase the tension on the brake cable, turn the screw counterclockwise; then - to keep this setting - turn the clamping button clockwise and hold down the adjustment screw.

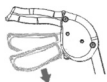
Fig.1 parking brake

Fig.2 braking

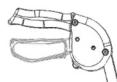
Fig.3 free movement

TIP: In case of any uncertainty, please contact Customer Service. Changing the brake settings can reduce the braking force.

MAINTENANCE

1. Please make sure ALL accessories are properly tightened at all times.

2. Please check that the screws, buttons and knobs of the wheels and moving wheels are tightened. Please check that the brake is functioning properly and that the handles are properly secured.

3. Worn or damaged parts must be replaced immediately.

CLEANING

To clean the walker, you need a damp cloth or sponge and warm soapy water. Normal disinfectants can also be used if required. Please make sure that the surfaces to be cleaned are not left with water or cleaning agents and that they are not slippery. Please do not use abrasive cleaning agents or scrapers as they may damage the paintwork.

Delivery range:

Walker, bag, instruction manual

METHOD OF DISPOSAL OF THE PRODUCT AFTER REMOVING FROM SERVICE

After the device is withdrawn from use, the medical device can be disposed of as normal household waste, except for electrical products - proceed in the manner appropriate for the disposal of electrical and electronic equipment.

This symbol indicates the user's maximum weight

SPECIFICATION:

width between handles 44 cm

seat height 63 cm

min. height of the handles 82 cm

maximum height of the handles 97 cm

total width 61 cm

total length 70 cm

seat size 22 cm

width when folded 22 cm

maximum load of the bag 5 kg

maximum user weight 136 kg

weight of a walker 9 kg

wheel diameter 29 cm, pneumatic wheels

ATTENTION: it is forbidden to use the product other than in accordance with its intended purpose

NOTE: Before use, check that all components are properly assembled.

ATTENTION: The manufacturer is not responsible for damage caused by neglect of maintenance, inadequate servicing, or resulting from failure to comply with the recommendations contained in this manual.

CAUTION: There may be a risk of tipping over if the product is used improperly. Please follow the boarding / dismounting / moving instructions.

CAUTION: When using and handling the product, as well as when assembling and adjusting the mechanisms, there may be a risk of entrapment and / or pinching of body parts of the user / accompanying person in openings / gaps between the elements. These steps should be performed with extreme caution. After completing the adjustment, stabilize the position by tightening the nuts / bolts carefully.

TARGET GROUP OF PATIENTS

Persons suffering from diseases, dysfunctions or injuries for the treatment/rehabilitation or compensation of which this device is intended (see section on intended use of the device in this instruction manual). The device can be purchased by the user on their own or on the recommendation of a doctor, therapist or other specialist. Whether you purchase the device yourself or on the recommendation of a doctor/therapist/other specialist, you must take into account the available sizes/necessary functions and variants of the device, the indications and contraindications for use and the information provided by the manufacturer.

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WARRANTY CARD

Model:
Serial number:

Point of sale stamp
and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the rollator, starting from the date of purchase by the Buyer. Products financed by the National Health Fund are covered with the extended warranty period, which covers half of the useful life specified in the Regulation on Guaranteed Services, provided that the useful life is specified in the Regulation in months or years.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, send it immediately to ANTAR (as specified in point 1) or contact the point of sale.
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum permissible weight of the user as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts subject to normal wear and tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED. IN ORDER TO EXTEND WARRANTY PERIOD FOR PRODUCTS FINANCED BY THE NATIONAL HEALTH FUND, ATTACH A COPY OF THE ORDER FOR THE SUPPLY OF MEDICAL DEVICES BEING ORTHOPAEDIC ITEMS AND AIDS.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net e-mail:antar@antar.net

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