AT51123

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EN

Thank you for purchasing our product, please read the instructions for use carefully. Do not use the product before reading and understanding this manual. If you do not understand the warnings, cautions and recommendations, please contact your health care professional or dealer to avoid product damage or personal injury.

CAUTION: Inspect all parts for damage caused in transit. If you notice such damage, DO NOT use the product. Contact the manufacturer for more information.

APPLICATION/INDICATIONS


The extremely comfortable and lightweight grabber for seniors is an excellent helper during ordinary activities. It is ideal as an aid for lifting objects from the floor as well as for taking things off high shelves. Its long and sturdy handle means you don't have to bend down or strain your hands. Uncomplicated and very easy to use, the grab bar for people with disabilities will increase the comfort of daily life. Limited range of motion, support for lifting objects from the floor or hard to reach places.

CONTRAINDICATIONS

Physical or mental limitations (such as visual impairment) that prevent safe handling.

- People who have limited sensation in their hands or arms may not adequately feel the weight or shape of an object grasped with the gripper, which may lead to dropping the object or using it incorrectly.
- People with impaired coordination may have difficulty controlling the gripper with precision, which may result in unintended movements or inability to grasp objects safely.
- If a person has limited strength in his or her hands or arms, he or she may not be able to use the gripper properly, especially if it requires more effort to squeeze or hold.
- Using a gripper, especially for people with balance limitations, may increase the risk of falls if the user must reach with it in an unstable position.
- For some neurological diseases, such as Parkinson's or multiple sclerosis, muscle tremors or spasms can make it difficult to use the gripper effectively and safely.

USE

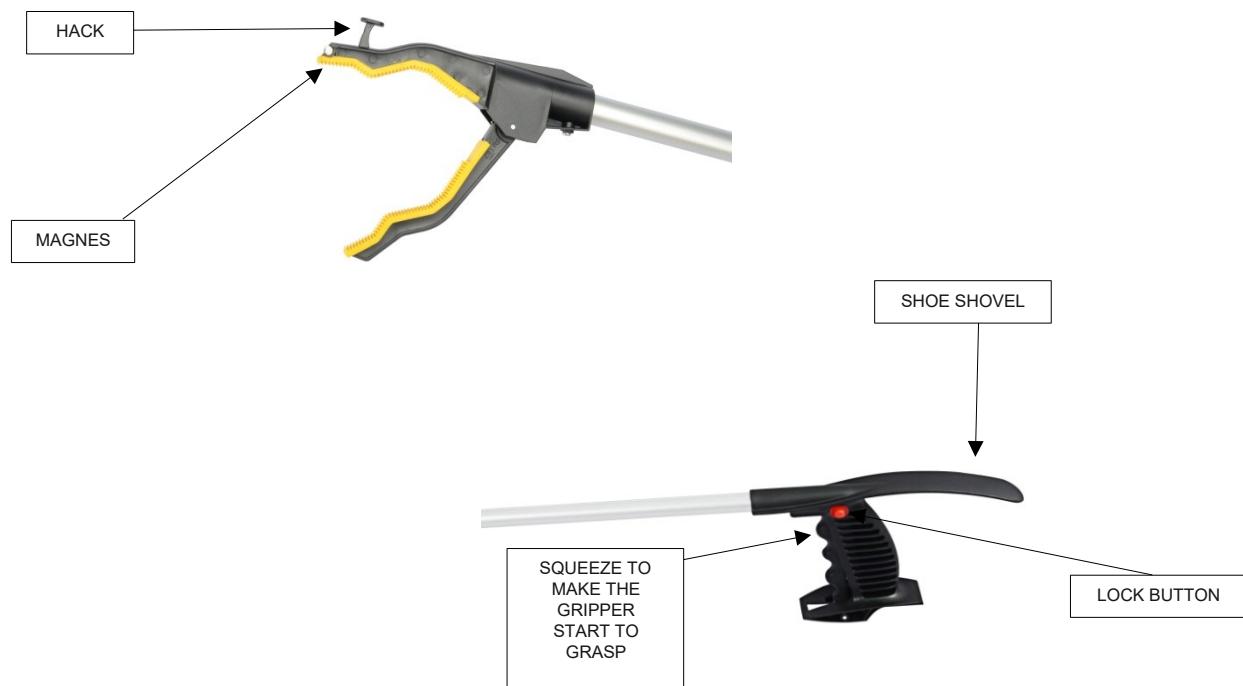FEATURES

- The top clamp head is equipped with magnets and a hook.
- Can attract metal objects and lift round objects.
- The handle is equipped with a flexible self-locking mechanism.
- The elongated upper end of the handle is in the shape of an arc, it assists when putting on shoes.

TECHNICAL SPECIFICATIONS

- Total length: 82.2cm
- Total width: 12.5cm
- Clip opening width: 9.3cm
- Total height: 4.5cm

- Total weight: 0.15kg
- Maximum object weight: 1kg

This mark indicates the maximum load capacity of the gripper

COMPONENTS:

Gripper, User manual

CAUTION:

In the event of a device-related “serious incident” that has directly or indirectly led to, may have led to, or may lead to any of the following:

- (a) death of a patient, user or other person, or
- (b) temporary or permanent deterioration of the health of a patient, user or other person, or
- (c) a serious threat to public health

the above “serious incident” must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

NOTES:

In case of pain, allergic reactions or other disturbing symptoms unclear to the user related to the use of the medical device, consult a health care professional.

NOTE: It is forbidden to use the product other than in accordance with its intended use.

NOTE: The manufacturer is not responsible for damage caused by negligent cleaning and maintenance, or as a result of failure to follow the recommendations in this manual.

TARGET PATIENT GROUP

Persons struggling with diseases, dysfunctions or injuries, for the treatment/rehabilitation or compensation of which this device is intended (see the section on the intended use of the device of this manual). The device may be purchased by the user independently or on the recommendation of a physician, therapist or other specialist. Whether purchasing the device on its own or as indicated by a doctor/therapist/other specialist, the available sizes/necessary functions and variants of the device, indications and contraindications for use, as well as information provided by the manufacturer, should be taken into account.

CLEANING AND CARE

The product should be cleaned by hand, using a cloth and water. A neutral detergent may be used. Do not clean with rough objects. Do not use strong detergents.

Remember to pay attention to the general condition of the product before each use. Does it show no signs of wear, no parts are missing, and all screws are properly tightened.

STORAGE

Store the product in a dry cool room, away from moisture and direct sunlight.

HOW TO DISPOSE OF THE PRODUCT

After the device is taken out of service, the medical device can be disposed of as ordinary municipal waste.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

 ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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