AT52201

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EN

Thank you for purchasing our product. Please read the user manual carefully before use. Do not use the product before reading and understanding this manual. If you do not understand the warnings, notes, or recommendations, contact a healthcare professional or the seller in order to avoid product damage or personal injury. **ATTENTION:** All parts should be checked for damage that may have occurred during transport. If such damage is noticed, **DO NOT** use the product. For more information, contact the seller.

TARGET PATIENT GROUP

Persons suffering from diseases, dysfunctions, or injuries for whose treatment, rehabilitation, or compensation this product is intended (see the section concerning the intended use of the product in this manual). The product may be purchased independently by the user or based on the recommendation of a doctor, therapist, or other specialist. In both cases—whether purchased independently or recommended by a healthcare professional—the available sizes, necessary functions, and product variants, as well as indications and contraindications for use and the information provided by the manufacturer, should be taken into consideration.

DESCRIPTION

The medical bed is intended for patient care in home environments or care facilities.

The structure of the bed is based on a steel frame divided into four sections, enabling adjustment of the patient's position. The bed base is made of welded steel, ensuring structural stability.

The headboard, footboard, and side rails are made of solid wood and MDF board. These elements provide protective and stabilizing functions during use.

The bed is equipped with an electric adjustment system controlled by a wired remote control. The adjustment allows changing the position of individual bed sections in order to adapt the patient's position to care or therapeutic needs.

The bed is also equipped with a lifting pole with a triangular handle, which assists the patient in repositioning in the bed.

For mobility, the bed is fitted with four braked castors.

The structure allows positioning in Trendelenburg and anti-Trendelenburg positions.

Full-length side rails can be raised or lowered depending on the user's needs.

The bed motor is equipped with a backup battery, which allows adjustments to be made in case of a power failure.

INTENDED USE

The electric rehabilitation bed is designed to assist elderly people, chronically ill patients, patients confined to a lying position, or those who have difficulty standing up or sitting down.

CONTRAINDICATIONS

Physical or mental limitations (e.g., visual impairment) that prevent the safe use of the product.

SAFETY

Do not assemble the product before reading this user manual.

Any modification of the product without the manufacturer's knowledge is the sole responsibility of the person performing the modification. Unauthorized modifications will void the warranty and invalidate the declaration of conformity with applicable directives.

It is forbidden to use the product in a manner other than described in this manual.

When the patient is not under supervision, the bed must be placed in the lowest position with the side rails raised to reduce the risk of falls and injuries.

Avoid improper use of the power cord. Hazards may arise from bending, cutting, or other mechanical damage.

Do not use the product for purposes other than those for which it was designed.

Keep this manual for reference or personnel training.

Warning! Do not exceed the maximum permissible load.

Warning! Always place the bed on a level surface.

Warning! Do not attempt to cross thresholds, steps, or other obstacles.

Warning! Do not place limbs near moving parts.

The bed must be equipped with appropriate elements to ensure that the power cord remains separated from moving parts or mechanisms of the medical bed during operation or movement in order to prevent damage to the cord.

Power cables must be adequately protected from damage caused by contact with moving parts or rubbing against sharp edges and surfaces inside the bed.

COMPONENTS: Hand control remote, Castors with individual braking system, Wooden headboard, Wooden footboard, Wooden side rails, Lifting pole (trapeze) with handle

TECHNICAL SPECIFICATIONS

Dimensions and adjustments

Overall length: 2147 mm

Overall width: 1028 mm

Overall height: 2273 mm

Mattress platform dimensions (bed frame): 1935 × 895 mm

Internal distance between headboard and footboard: 2045 mm

Height adjustment range (floor to bed frame)

Highest position: 828 mm

Lowest position: 418 mm

Section adjustment range

Maximum backrest inclination: 70°

Maximum leg section inclination: 34°

Maximum foot section inclination: 22°

Trendelenburg position: ±12°

Reverse Trendelenburg position: ±12°

Mattress platform sections

Head section: 800 × 740 mm
 Fixed section: 895 × 206 mm
 Leg section: 800 × 300 mm
 Foot section: 800 × 465 mm
 Maximum height of leg section from bed frame: 230 mm
 Mobility system
 Castors: 4 pcs, diameter 4 inches (approx. 102 mm)
 Type: swivel with brake
 Electrical components
 Remote control: 1 pc
 Electric motors: 2 pcs (dual motor)
 Control box: 1 pc
 Power supply: 100–240 V
 Duty cycle: 2 minutes operation / 18 minutes rest
 Operating parameters
 Maximum user weight: 135 kg
 Maximum safe load (patient + mattress + accessories): 170 kg
 Maximum safe working load of lifting mechanism: 270 kg
 Maximum mattress weight: 20 kg
 Maximum accessories weight: 25 kg
 Bed weight: 95 kg
 Operating conditions
 Operating temperature: -10°C to +50°C
 Storage temperature: -10°C to +50°C
 Operating environment: indoor use
 Materials
 Frame: steel
 Headboard / footboard: solid wood + MDF
 Side rails: solid wood
 Castors: steel

This symbol indicates the maximum user weight.

ATTENTION!

Individual product dimensions may vary within ±2 cm.

SWIVEL CASTORS WITH BRAKE

The castors are intended only for moving the bed when necessary. At all other times they must remain locked.

To lock the wheels, press the brake downward.

To unlock the wheels, lift the brake to the horizontal position.

ATTENTION!

During bed assembly, the wheels at the headboard and footboard must always remain locked. Lock them by pressing the brake pedal with your foot.

BED BASE ASSEMBLY

Remove the headboard and the frame module of the orthopedic bed intended for the backrest function from the packaging.

Insert the frame into the bayonet mounts located at the headboard.

Pull the spring lock button towards the inside of the bed and turn it. Then hook the frame into the bayonet mounts. Finally, pull the button towards the inside of the bed again and turn it until it locks.

Remove the footboard and the bed frame module intended for the pelvis and legs section from the packaging, and repeat the same steps performed during the headboard assembly.

Join the two previously assembled bed base modules together and connect them in the central section on both sides to form the base of the orthopedic bed.


Remove the four black clips from the packaging and screw them into the appropriate holes in the central part of the base on both sides to ensure the stability and safety of the assembled bed structure.

MOTOR ASSEMBLY

Remove the bed motor from the packaging and insert the two backup batteries into the appropriate slots.

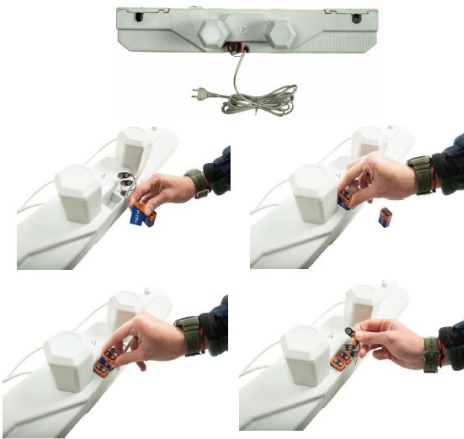Insert the blanking plugs into both batteries, ensuring they are correctly placed in the appropriate slots.


The motor has "Head" and "Feet" markings on both ends to facilitate correct orientation: the "Head" marking should face the headboard, the "Feet" marking should face the footboard.


Attach the motor to the frame by pushing it from the bottom upward until the motor locks correctly onto the frame.

Connect the headboard and footboard motor pins to their designated sockets.

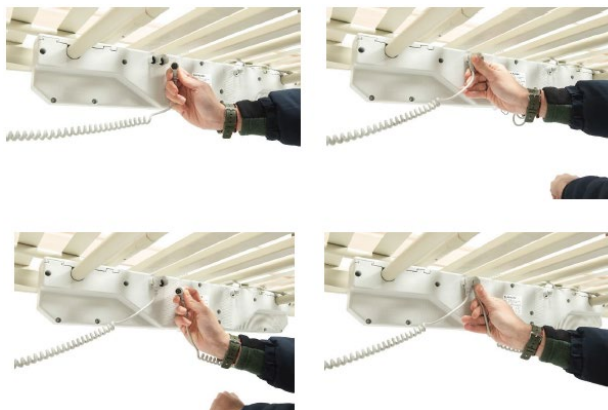Plug the remote control connector into the remote control socket.

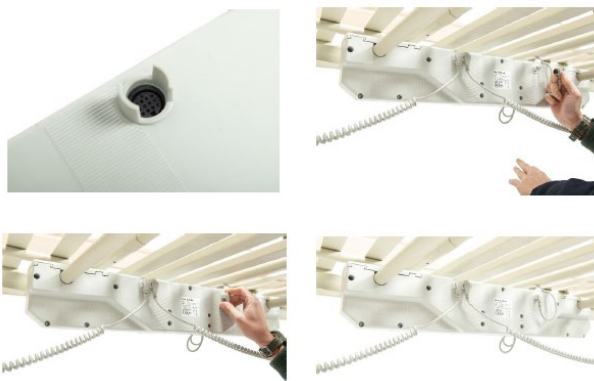Remove the two locking elements from the packaging. Insert them into the two upper openings of the motor until they are fully inserted and locked.

SIDE RAIL ASSEMBLY

Remove the side rails and guides from the packaging.
At the ends of the rails, insert two guides into the appropriate holes.

ATTENTION!

Raise the bed base before starting the installation of the side rails.
Press the grey button and insert the rails with the installed guides from bottom to top until you hear a click indicating they are locked in the final position.

The side rails are raised by pressing the grey button and guiding them upward until they engage in the safe position. To lower them, press the button again and guide the rails downward in a controlled manner.

ASSEMBLY OF THE LIFTING POLE AND TRIANGLE HANDLE

Remove the lifting pole and the triangular handle from the packaging. Insert the hook of the triangular handle into the hole located at the top of the lifting pole.

Insert the lifting pole into the designated socket in the bed, making sure it is securely and stably mounted.

REMOTE CONTROL FUNCTIONS

1. Raising / lowering the backrest
2. Raising / lowering the leg section
3. Simultaneous raising / lowering of the leg and foot sections
4. Raising / lowering the bed frame
5. Trendelenburg / anti-Trendelenburg position

IMPORTANT

On the back side of the remote control there is a key lock system. The remote control can be locked or unlocked using the key integrated in the remote control cable.

To lock the buttons:

Turn the key until the arrow aligns with the closed padlock symbol.

To unlock the buttons:

Turn the key until the arrow aligns with the open padlock symbol.

Problem — Solution

The bed does not respond to up/down commands

Check whether the plug is properly connected

Check and clean the contacts

Reduce the load

Remote control LED does not work

Check whether the connector is properly connected

The bed moves irregularly

Clean the inside of the remote control

Contact the manufacturer / retailer

Lubricate the joints

Headboard / footboard are not properly installed

Contact the manufacturer / retailer

Metallic sounds during movement

Contact the manufacturer / retailer

ATTENTION

If a serious incident related to the device occurs that directly or indirectly led, might have led, or may lead to any of the following:

a) death of a patient, user, or another person

b) temporary or permanent serious deterioration of the health of a patient, user, or another person

c) a serious public health threat

the incident must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides.

ATTENTION: If pain, allergic reactions, or other disturbing or unclear symptoms related to the use of the medical device occur, consult a healthcare professional.

ATTENTION: It is prohibited to use the product in a manner other than in accordance with its intended purpose.

ATTENTION: During the use, handling, assembly, or adjustment of the mechanisms of the product there is a risk of trapping and/or crushing parts of the body of the user or caregiver in openings or gaps between the elements. These activities must be carried out with particular caution. After adjustment, stabilize the position by tightening nuts and screws securely.

CLEANING AND MAINTENANCE

Before each use of the bed, check the wheel brakes, the rail fixings, and the boards at the head and foot of the bed.

The bed should be stored in a dry and ventilated room.

Keep the bed clean. Clean using a damp cloth.

Before moving the bed, place it in the lowest position. Always use the headboard or footboard to move the bed. Do not use the side rails to move the bed. Exercise caution when maneuvering the bed. After moving the bed to the desired location, lock the wheels.

STORAGE: The product should be stored in a dry place and should not be exposed to direct sunlight.

DISPOSAL: After the product is withdrawn from use, it should be disposed of as regular municipal waste. Remember to follow waste segregation regulations. Protect the environment.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp. z o.o. ,03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.
10. The warranty does not cover parts and components subject to natural wear and tear, such as wheels, tires, rubber pads, fuses, etc., the wear of which results from normal use of the product. The warranty period for batteries is 6 months.

ANTAR MEDICAL Sp. z o.o.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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