
AT03104/AT03015/AT03106
INSTRUKCJA UŻYTKOWNIA, MATERAC PRZECIWODLEŻYNOWY PIANKOWY
USER MANUAL, FOAM ANTI-BEDSORE MATTRESS
BEDIENUNGSANLEITUNG, SCHAUMSTOFF-ANTI-DEKUBITUS-
MATRATZE
NÁVOD K POUŽITÍ, PĚNOVÁ MATRACE PROTI PROLEŽENINÁM
NÁVOD NA POUŽITIE, PENOVÝ PROTIVREDOVÝ MATRAC
MANUEL D'UTILISATION, MATELAS ANTI-ESCARRES EN MOUSSE
GEBRUIKSAANWIJZING, SCHUIM ANTI-DECUBITUS MATRAS
MANUAL DE INSTRUCCIONES, COLCHÓN ANTIÚLCERAS DE ESPUMA
MANUALE DI ISTRUZIONI, MATERASSO ANTI-DECUBITO IN SCHIUMA
BRUKSANVISNING, SKUM ANTI-LIGGSÅRSMADRASS

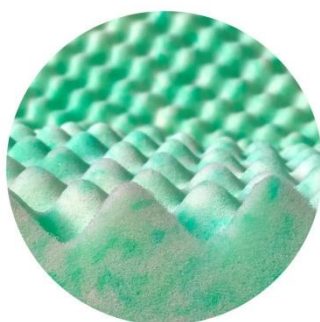PL

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EN

Thank you for purchasing our product, please read the instructions for use carefully. Do not use the product before reading and understanding this manual. If you do not understand the warnings, cautions and recommendations, please contact your health care professional or dealer to avoid product damage or personal injury.

CAUTION: Inspect all parts for damage caused in transit. If you notice such damage, DO NOT use the product. See your dealer for more information.

TARGET PATIENT GROUP

Persons struggling with diseases, dysfunctions or injuries, for the treatment/rehabilitation or compensation of which this product is intended (see the section on the purpose of the product of this manual). The product may be purchased by the user independently or on the recommendation of a doctor, therapist or other specialist. Whether purchasing the product on its own or as indicated by a doctor/therapist/other specialist, the available sizes/necessary functions and variants of the device, indications and contraindications for use, as well as information provided by the manufacturer, should be taken into account.

APPLICATION/INDICATIONS

Anti-decubitus foam mattress AT03104/AT03105 is a modern solution designed for people who spend most of their time lying down due to limited mobility. Thanks to the advanced double-layer construction, this mattress not only provides comfort, but also effectively prevents the formation of bedsores.

For people with limited mobility: The mattress meets the needs of patients who spend most of their time in bed, both at home and in medical centers.

Comfort and prevention of bedsores: Thanks to special foam, the mattress reduces pressure on the body while providing comfort.

Ergonomic fit: The design supports correct posture when lying down, which benefits the health of the spine.

Durability: The mattress is designed to accommodate users weighing up to 130 kg, making it suitable for a wide range of patients.

Two layers of foam:

The first layer of polyurethane foam (10 cm) offers firm support and resilience.

The second layer, made of memory effect foam (2 cm), adapts perfectly to the body, which reduces pressure on tissues.

Choice of customized covers:

The waterproof cover protects the mattress from moisture, which increases its durability and makes it easier to keep clean.

The disinfectable cover is easy to clean, which is crucial in a medical environment.

The AT03106 anti-decubitus foam mattress is an indispensable solution for people with limited mobility who spend most of their time lying down. Thanks to innovative technologies and carefully selected materials, it not only provides effective protection against bedsores, but also optimal comfort.

Advantages of using mattress AT03106

Effective prevention of bedsores: the mattress is ideal for the care of people lying down, preventing the formation of painful skin lesions.

Convenience of use: Thanks to its ergonomic design and anti-allergenic properties, it provides complete comfort for both patients and their caregivers.

Durability and hygiene: Waterproof cover and deformation-resistant foam guarantee long-term use.

Versatility of use: The mattress is suitable for people weighing up to 100 kg, making it a versatile choice for a variety of patients.

Thanks to its modern design and the materials used, the AT03106 anti-decubitus foam mattress not only helps prevent pressure sores, but also promotes healthy sleep. Moreover, its ergonomic structure adapts to the shape of the body, offering maximum comfort to the user.

FEATURES AND TECHNICAL SPECIFICATIONS OF MATTRESS AT03104/AT03105:

Mattress AT03104/AT03105:

Dual-layer foam for greater comfort and support:

Base layer: 10 cm of polyurethane foam, which is responsible for stability and proper body support while sleeping.

Top layer: 2 cm of memory effect foam, which allows the mattress to conform to the shape of the body to reduce pressure on sensitive areas of the skin.

Optimal body weight distribution:

The mattress is designed to distribute the user's weight evenly, which helps maintain the natural curves of the spine. This not only improves comfort, but also reduces the risk of pain.

Various covers:

Waterproof: Effectively protects the mattress from moisture, which translates into durability and hygiene.

Disinfectable: Easy to clean, ideal for use in medical facilities where cleanliness is a priority.

Dimensions AT03104: 200 cm x 90 cm x 12 cm

Dimensions AT03105: 200 cm x 120 cm x 12 cm

Maximum user weight: 130 kg

This symbol indicates the maximum weight of the user!

FEATURES AND TECHNICAL SPECIFICATIONS AT03106

Anti-allergenic polyurethane foam with castor oil

The mattress is made of polyurethane foam, which is safe for allergy sufferers. Moreover, the addition of castor oil has an antibacterial effect, ensuring hygienic conditions of use.

Special piping structure

The unique piped surface of the mattress allows for better air circulation. This reduces the risk of dampness, which effectively protects against the formation of bedsores.

Even distribution of body weight

Thanks to its structure, the mattress evenly distributes body weight, which helps prevent pressure points. As a result, it supports the natural curves of the spine and minimizes the risk of skin damage.

Waterproof and easy-to-clean cover

The waterproof cover protects the mattress from moisture and dirt. Moreover, it is easy to clean, which increases the hygiene and durability of the product.

Dimensions AT03106: 200 cm x 90 cm x 5 cm

Material: polyurethane foam with castor oil

Maximum user weight: 100 kg

This symbol indicates the maximum weight of the user!

CONTRAINDICATIONS

Physical or mental limitations (e.g., visual impairment) that prevent safe handling of the product.

Significant body deformities - Individuals with severe anatomical deformities may not receive adequate support, which may lead to discomfort or impairment.

Allergies to materials used in the mattress - If you are allergic to latex, polyurethane foam or other ingredients used in the mattress, its use may cause allergic reactions.

Severe spinal or pelvic injuries - Foam mattresses may not provide sufficient stability for patients with severe orthopedic or neurological injuries.

Advanced stage of decubitus ulcers - In patients with pre-existing, deep decubitus wounds, a foam mattress may not provide adequate pressure relief and a more advanced decubitus system (such as an air mattress) may be required.

Circulatory failure - Patients with circulatory disorders require specialized monitoring and variable pressure mattresses that can improve blood circulation.

Inability to change body position - In patients who are unable to change body position on their own and require frequent repositioning, a foam mattress may not be sufficient.

Overweight or obesity - A foam mattress may not provide adequate support for patients whose weight exceeds its maximum load, which may lead to deformation or damage.

LIST OF COMPONENTS

Mattress with cover, instruction manual

CAUTION:

In the event of a device-related “serious incident” that directly or indirectly led to, may have led to, or may lead to any of the following:

- (a) death of a patient, user or other person, or
- (b) temporary or permanent deterioration of the health of a patient, user or other person, or
- (c) a serious threat to public health

the above “serious incident” must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

NOTES:

In case of pain, allergic reactions or other disturbing symptoms unclear to the user related to the use of the medical device, consult a health care professional.

CAUTION: It is forbidden to use the product other than in accordance with its intended use.

CLEANING

Washing instructions for the waterproof cover (standard):

Hand wash the cover only.

Clean using a damp sponge with laundry detergent.

Maximum washing temperature 40°C.

Do not bleach! (do not use free chlorine emitting compounds for bleaching).

Do not dry clean!

Do not tumble dry.

Do not iron! (do not allow contact with heated surfaces such as a radiator).

Deterioration of visual-aesthetic properties (type of surface wrinkling) does not affect change the functional parameters of the material (physical-mechanical properties).

Washing instructions for disinfecting cover:

Do not wash the foam!

Clean using a damp sponge with detergent.

Maximum washing temperature 95°C.

Do not bleach! (do not use free chlorine emitting compounds for bleaching).

Do not dry clean!

Do not tumble dry.

Do not iron! (do not allow contact with heated surfaces such as a radiator).

Deterioration of visual-aesthetic properties (type of surface wrinkling) does not affect change the functional parameters of the material (physical-mechanical properties).

Cleaning foam requires a delicate approach to avoid damaging its structure. Here is a guide on how to clean it properly:

What you will need:

Mild detergent or dishwashing liquid

Warm water

Baking soda (optional for deodorizing)

Vacuum cleaner (preferably with a soft brush)

Spray bottle

Soft cloth or sponge

Enzymatic agent (optional for stain removal)

Cleaning steps:

Vacuum the foam:

Remove dust and debris: Use a vacuum cleaner with a soft brush to gently vacuum the foam surface. This will remove dust, dirt and debris before wet cleaning.

Spot clean with a mild detergent:

Mix detergent: Combine a few drops of mild detergent with warm water in a spray bottle.

Gently spray the foam: Avoid soaking the foam. Lightly spray the mixture on the soiled area.

Dry the stains: Use a soft cloth or sponge to gently dry the area. Do not rub, as this may damage the structure of the foam.

Rinse the area: If necessary, spray lightly with clean water to remove any residual detergent.

Dealing with difficult stains:

Biological stains or urine: Use an enzymatic agent. Lightly apply the product to the stain, let it sit for a while (according to the product's instructions), and then dry with a clean cloth.

Greasy stains: Use a small amount of dishwashing liquid or a special upholstery cleaner designed for foam materials.

Deodorization (optional):

Use baking soda: Sprinkle a layer of baking soda on the foam and let it absorb odors for several hours or even overnight.

Dust off the baking soda after enough time has passed.

Drying the foam:

Completely air-dry: Foam is very absorbent, so it may take some time to dry. Leave it in a well-ventilated area.

Use a fan to speed up the process: Never use a direct heat source, such as a hair dryer, as this can damage the foam. It is best to use a cool fan.

Make sure the foam is completely dry before reusing it to prevent mold or mildew growth.

What to avoid:

Never immerse the foam in water, as it can be difficult to dry and may be damaged.

Avoid using strong chemicals such as bleach or alcohol-based cleaners, which can weaken the structure of the foam.

No direct heat: Heat can damage the structure of the foam, so avoid drying in a dryer or using hot air.

STORAGE

Store the product in a dry place, do not expose to direct sunlight.

UTILIZATION

After decommissioning, treat the product as ordinary municipal waste. Be sure to comply with waste segregation standards. Take care of the environment.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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