## **AT53085**

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**ORTEZA KOLANA**

**INSTRUCTIONS FOR USE**

**KNEE ORTHOSIS**

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**KNIE-ORTHESE**

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**NÁVOD NA POUŽITIE**

**ORTÉZA KOLENA**

**MODE D'EMPLOI**

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**KNEE ORTHOSIS**

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**ORTOSIS DE RODILLA**

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**ORTESI DEL GINOCCHIO**

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**KNÄ-ORTOS**

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En

Thank you for purchasing our product, please read the instructions for use carefully. Do not use the product before reading and understanding this manual. If you do not understand the warnings, cautions and recommendations, please contact your health care professional or dealer to avoid damage to the product and bodily injury.

**NOTE:** Inspect all parts for damage caused in transit. If you notice such damage, DO NOT use the product. Contact the manufacturer for more information.

#### USE/INDICATIONS

- acute/chronic pain syndromes
- treatment of knee joint injuries
- Injuries and dysfunctions without immobilizing the joint
- support of the ligamentous apparatus of the knee joint
- post-operative care of cartilage or meniscus

#### CONTRAINDICATIONS

Physical or mental limitations (e.g., visual impairment) that prevent safe handling.

Any injury or illness that may affect the comfort and safety of using the product.

-skin infections: If there are infections, ulcers or other skin conditions on the skin around the knee, the use of the orthotic may worsen the skin condition or cause further infections.

-lack of adequate tolerance: Some patients may experience discomfort or pain when wearing the orthosis, which may affect their ability to engage in normal physical activity.


-muscle atrophy: In people with advanced muscle atrophy in the knee area, the use of an orthotic may be ineffective or lead to further muscle weakness.

-Serious joint injuries: For severe joint injuries, such as fractures or ligament ruptures, the use of an orthotic may require specialized evaluation and may not be sufficient for treatment.

-unstable joint: If the joint is very unstable or there is a risk of further damage, the use of an orthotic may not be sufficient or surgical intervention may be required.

#### SIZE

| SIZE | THIGH (CM) | CALF (CM) |
|------|------------|-----------|
| S    | 39-47      | 29-37     |
| M    | 47-54      | 37-43     |
| L    | 54-60      | 43-50     |
| XL   | 60-67      | 50-57     |

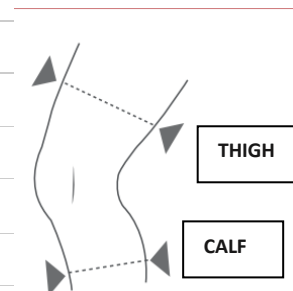**\*Product comes in right and left leg versions.**

#### FEATURES

- Full structure provides good stability
- soft lining
- The orthosis can be used at home or for sports.
- product anti-slip
- product easy to take off and put on
- ROM hinges with plug-in bend/extension stops.
- Orthosis of the knee type OA

According to the ergonomic design of the knee joint, the curve of the bracket adapts to the shape of the leg, reducing pressure.

There are metal covers on the front and back of the product, and brackets are used to support the left and right sides to increase protection and ensure a proper fit.

High-quality Velcro straps are durable and easy to fasten. Suitable for patients of different body types. The inner side of the strap is attached with a non-slip silicone strip, which not only improves the fit of the wearer, but also prevents the strap from falling off easily during exercise. At the same time, it makes the wearing experience softer and more comfortable.

**COMPONENT PARTS:** Knee brace, Instruction manual, hinges, allen wrench

## NOTES:

In the event of a product-related "serious incident" that directly or indirectly led to, may have led to, or may lead to any of the following:

- (a) death of a patient, user or other person or
- (b) temporary or permanent deterioration of the health of the patient, user or other person or
- (c) a serious threat to public health

the above "serious incident" should be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

## NOTES:

If you experience pain, allergic reactions, or other disturbing symptoms that are unclear to you related to the use of the medical device, consult a healthcare professional.

**NOTE: it is forbidden to use the product other than for its intended purpose .**

**NOTE:** The manufacturer is not responsible for damage caused by negligent cleaning and maintenance, or resulting from failure to follow the recommendations in this manual.

## TARGET PATIENT GROUP

Persons struggling with diseases, dysfunctions or injuries, for the treatment/rehabilitation or compensation of which this device is intended (see the section on the purpose of the device of this manual). The device may be purchased by the user independently or on the recommendation of a physician, therapist or other specialist. Whether purchasing the product on its own or as indicated by a doctor/therapist/other specialist, the available sizes/necessary functions and variants of the device, indications and contraindications for use, as well as information provided by the manufacturer, should be taken into account.

## USE

1. unfasten all straps(velcro).
2. place the axis of the hinge in line with the axis of rotation of the knee. (FIG 1)
- Wear a knee brace when the patient is sitting and the knee joint flexion is 80-90 degrees.**
3. press the stabilizer against the leg and fasten the first strap behind the calf. Adjust the length of the strap to the desired fit using the Velcro fastener. Then fasten the strap around the lower calf muscles. (FIG 2)
- 4 Attach the remaining straps in the same manner as the previous straps (FIG 3)
5. stand and check the tension of the straps (FIG 4)
6. Walk around with the orthosis on to see if the strap causes discomfort or exacerbates the pain. If so, readjust the strap in a sitting position, then stand and walk again to see if the pain has subsided.
7. adjust the hinge with an allen wrench (FIG 5)

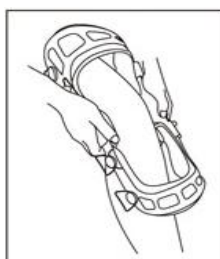
Fig 1

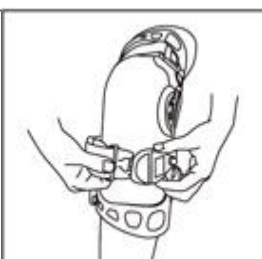
Fig 2

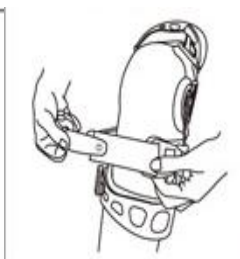
Fig 3

Fig 4

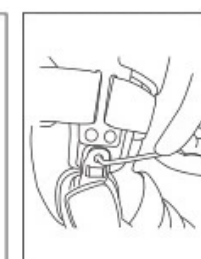
Fig 5

## HOW TO ADJUST THE HINGE

1. use an allen wrench to remove the screws and remove the hinge cover.
2. Remove the existing stop and insert new plug stops that have the desired degree of flexion/extension.
3. Put the hinge cover and screws back on. Make sure that the screws are tightened well. Both the medial (inner) and lateral (outer) hinge stops must be identical in degree.

Limitation of uprightness at: 0°, 10°, 20°, 30°, 40°.



Limitation of bending at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°.

Limitation of immobilization at: 0°, 10°, 20°, 30°.

## CLEANING AND CARE

**The fabric parts of the orthotic can be detached and should be hand washed** using cold water in sunlight, hand wash using neutral detergent. Note: Do not wash the metal parts of the orthosis!

- Do not iron.
- Do not dry-clean.
- Do not use bleach.
- Protect from moisture.
- Do not expose to high temperatures.
- Do not tumble dry.
- Dry unfolded, gently squeeze before drying.
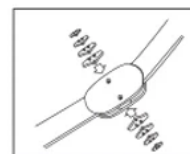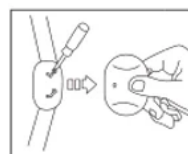
|  |                                       |
|--|---------------------------------------|
|  | HAND WASH AT UP TO 40 DEGREES CELSIUS |
|  | Do not use bleach                     |
|  | Do not iron                           |
|  | Do not tumble dry                     |
|  | Do not dry-clean                      |

Remember to regularly check that the screws are well screwed in and that there is no play, which can affect the safe and comfortable use of the product.

Do not use the product when it is contaminated.

#### STORAGE

Store the product in a dry cool room, away from moisture and direct sunlight.

#### METHOD OF DISPOSAL OF THE PRODUCT

Once the device is decommissioned, the medical device can be disposed of as ordinary municipal waste.

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Model: .....  
Serial number: .....

Point of sale stamp

and a legible signature of the Seller: .....

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.  
03-068 Warsaw, ul. Zawiślańska 43  
Tel. 22 518 36 00, Fax 22 518 36 30  
www.antar.net

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