**INSTRUKCJA UŻYTKOWANIA AT53034 ORTEZA STAWU
KOLANOWEGO Z REGULACJĄ ZAKRESU RUCHOMOŚCI
INSTRUCTION MANUAL AT53034 KNEE JOINT ORTHOSIS WITH RANGE
OF MOTION ADJUSTMENT
GEBRAUCHSANWEISUNG AT53034 KNIESORTHESE MIT
EINSTELLBAREM BEWEGUNGSUMFANG
NÁVOD K POUŽITÍ AT53034 KOLENNÍ ORTÉZA S NASTAVITELNÝM
ROZSAHEM POHYBU
NÁVOD NA POUŽITIE AT53034 KOLENNÁ ORTÉZA S NASTAVITEĽNÝM
ROZSAHOM POHYBU
MODE D'EMPLOI AT53034 GENOUILLÈRE À AMPLITUDE DE
MOUVEMENT RÉGLABLE
GEBRUIKSAANWIJZING AT53034 KNIEBRACE MET INSTELBAAR
BEWEGINGSBEREIK
INSTRUCCIONES DE USO AT53034 RODILLERA CON AMPLITUD DE
MOVIMIENTO AJUSTABLE
ISTRUZIONI PER L'USO AT53034 GINOCCHIERA CON RANGE DI
MOVIMENTO REGOLABILE
BRUKSANVISNING AT53034 KNÄSTÖD MED JUSTERBART
RÖRELSEOMFÅNG**

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En

Thank you for purchasing our product, please read the instructions for use carefully. Do not use the product before reading and understanding this manual. If you do not understand the warnings, cautions and recommendations, please contact your health care professional or dealer to avoid damage to the product.

NOTE: Inspect all parts for damage caused in transit. If you notice such damage, DO NOT use the product. Contact the manufacturer for more information.

WARNINGS AND CAUTIONS

- 1 This device is intended for one patient only. It is not intended for reuse
By another patient.
- (2) If additional pain, swelling or other symptoms occur while using the device, discontinue
Use and contact your doctor.
3. this device is designed to protect the injured area within a limited range of motion in accordance with the
with a rehabilitation program or other rehabilitation measures prescribed by a doctor.
- 4 Storage: in a dry and clean environment, well ventilated, without corrosive gases,
At room temperature.
- 5 Regularly check the equipment for damage and wear. In case of any damage
Stop using the product immediately.
6. repair and disinfection should be performed only by authorized persons!
7. check that all screws are properly tightened before each installation.
8. during use and operation and during its assembly and adjustment of mechanisms may occur
Danger of entrapment and/or compression of body parts of the user/companion in the openings/
The gaps between the elements.

USE/INDICATIONS

The product is aimed at people who, as a result of diseases or various ailments and injuries, need to stabilize the knee joint.

- reconstruction of ligaments: ACL, PCL
- disorders of the meniscus
- stability of the knee joint
- prompt
- post-traumatic rehabilitation process (e.g., after sprains, patellofemoral conflict)
- limited range of motion after reconstructive surgery

CONTRAINDICATIONS

Physical or mental limitations (e.g., visual impairment) that prevent safe handling.

FEATURES

- possibility to adjust the bending angle as desired (0°-90°)
Restriction of extension at: 0°, 10°, 20°, 30°, 40°.
- Restriction of flexion at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°.
- Hypoallergenic silicone strips that provide a soft feeling on the skin, while protecting the
camera from slipping
- Additional knee stabilizing inserts on both sides
- Breathable padding combines good cushioning and optimal heat transfer.
- color: black

SIZE

NOTE: Orthosis comes in right and left leg version

In order to fit the right size, you should:

Measure the abdomen of the thigh 15 cm above the knee:

S: 39 - 47 cm, M: 47-53 cm, L: 53-60 cm, XL: 60-67 cm

Measure the calf abdomen 13 cm below the knee:

S: 29 - 37 cm, M: 37 - 43 cm, L: 43 - 50 cm, XL: 50 - 57 cm

SCOPE OF DELIVERY:

Orthosis, Instructions for use

NOTES:

In the event of a product-related "serious incident" that directly or indirectly led to, may have led to, or may lead to any of the following:

- (a) death of a patient, user or other person or
- (b) temporary or permanent deterioration of the health of the patient, user or other person or
- (c) a serious threat to public health

the above "serious incident" should be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

NOTES:

If you experience pain, allergic reactions, or other disturbing symptoms that are unclear to you related to the use of the medical device, consult a healthcare professional.

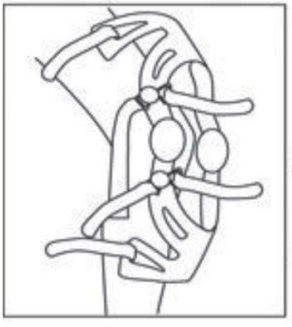
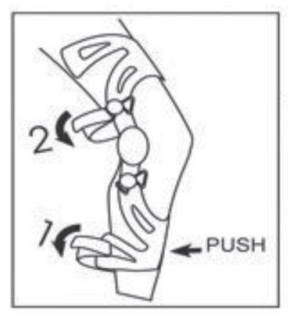
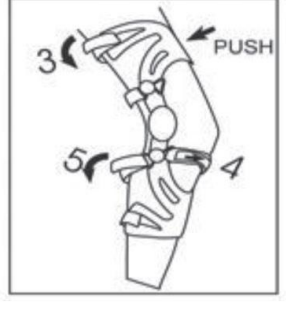
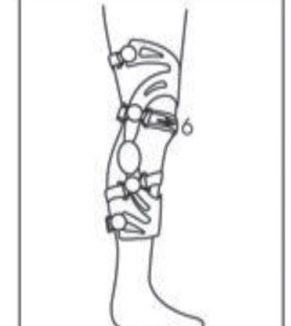
NOTE: it is forbidden to use the product other than for its intended purpose .

NOTE: The manufacturer is not responsible for damage caused by negligent cleaning and maintenance, or resulting from failure to follow the recommendations in this manual.

TARGET PATIENT GROUP

Persons struggling with diseases, dysfunctions or injuries, for the treatment/rehabilitation or compensation of which this device is intended (see the section on the purpose of the device of this manual). The stabilizer may be purchased by the user independently or on the recommendation of a physician, therapist or other specialist. Whether purchasing the orthosis on its own or as indicated by a doctor/therapist/other specialist, the available sizes/necessary functions and variants of the device, indications and contraindications for use, as well as information provided by the manufacturer, should be taken into account.

USE

1 Before applying, measure the circumference of the thigh and calf. Using the chart, choose the correct size of the orthotic. 2. undo all fasteners. The stabilizer should be worn in a sitting position. Bend the knees at an angle of about 30 degrees. 3. place the stabilizer on the knee joint. The center hinge should be placed in line with the knee joint.		4. rest the lower part of the stabilizer against the shin and tighten (1). Then fasten the attachment above the knee (2).	
5. make sure that the outer part of the stabilizer is firmly attached to the leg, then fasten the attachment (3) on the back of the thigh. Then, fasten the front (4) and rear (5) fasteners under the lower part of the knee. Be sure to fasten the fasteners in the correct order.		6. In the standing position, fasten the fastener above the knee (6). Note: Be sure to fasten the mounts in the order described in the instructions. Tighten each fastening strap so that the fastening is stable but does not restrict movement. After 10-15 minutes of activity, all fastenings should be loosened and then re-tighten in the correct order. This is necessary to limit swelling and excessive pressure on tissues.	

NOTE: a Allen key is required to adjust the bending angle

CLEANING AND CARE

The orthosis can be cleaned by hand with a damp washcloth and mild detergent. Avoid corrosive substances. Do not use rough objects for cleaning.

The orthosis has a sponge material that can be removed and washed by hand using mild soap. Once the material is dry, it can be attached back to the orthosis. Do not tumble dry, use bleach and do not dry clean.

STORAGE

Store the product in a dry cool room, away from moisture and direct sunlight.

METHOD OF DISPOSAL OF THE PRODUCT

Once the device is decommissioned, the medical device can be disposed of like regular municipal waste

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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