AT03023

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The cushion was developed for carers who take care of people, who are in bed permanently or for long periods of time. By using the cushion, turning the patient will not require a lot of force and for the person lying down this activity will be much more comfortable.

FEATURES:


- Resistant to deformation
- Can be surface disinfected
- Dimensions: 49 cm x 41.5 cm x 10.6 cm

USE/INDICATIONS

- elderly, disabled or paralysed people
- obese patients
- during general anaesthesia and other cognitive dysfunctions
- physical dysfunctions

USE

The cushion should be slid under the patient's legs and then turned in the desired direction (you can keep the patient in a particular position by attaching the cushion to the bed side with an additional attachment strap).

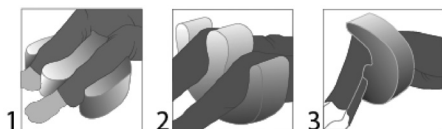CLEANING AND MAINTENANCE

The cushion should be cleaned using a damp cloth with detergent. Machine washable up to 95°C. Do not tumble dry, iron, bleach or dry clean. Dry the unfolded product in the shade - do not wash the foam. The cover can be surface disinfected and then wiped with a clean cloth.

CONTRAINDICATIONS

Physical or mental limitations (e.g. visual impairment) that prevent safe handling.

STORAGE

Store the product in a dry cool room, away from moisture and direct sunlight.

TARGET PATIENT GROUP

Persons suffering from diseases, dysfunctions or injuries for the treatment/rehabilitation or compensation of which this device is intended (see section on intended use of the device in these instructions). The device can be purchased by the user on their own or on the recommendation of a doctor, therapist or other specialist. Whether you purchase the device yourself or on the recommendation of a doctor/

therapist/other specialist, you must take into account the available sizes/necessary functions/size/variants of the device, the indications and contraindications for use, as well as the information provided by the manufacturer.

NOTE: It is prohibited to use the product other than for its intended purpose.

ATTENTION: The manufacturer is not responsible for damage caused by negligent maintenance, inadequate servicing or resulting from failure to observe the recommendations contained in this manual.

CAUTION:

In the event of a device-related „serious incident“ that directly or indirectly led to, could have led to or may lead to any of the following:

- (a) death of a patient, user or other person, or
- (b) temporary or permanent deterioration of the health of a patient, user or other person, or
- (c) a serious risk to public health

the above „serious incident“ must be notified to the manufacturer and to the competent authority of the Member State in which the user or patient is resident. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

NOTES:

In the event of pain, allergic reactions or other disturbing symptoms unclear to the user related to the use of the medical device, consult a healthcare professional.

HOW TO DISPOSE OF THE DEVICE AT END-OF-LIFE

Once the device has been withdrawn from service, the medical device can be disposed of as normal municipal waste.

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WARRANTY CARD

Model:
Serial number:

Point of sale stamp
and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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