AT53073
INSTRUKCJA UŻYTKOWNIA, PROSTOTRZYMACZ
USER MANUAL, POSTURE CORRECTOR
BEDIENUNGSANLEITUNG, HALTUNGSKORREKTUR
NÁVOD K POUŽITÍ, KOREKTOR DRŽENÍ TĚLA
NÁVOD NA POUŽITIE, KOREKTOR DRŽANIA TELA
MODE D'EMPLOI, CORRECTEUR DE POSTURE
GEBRUIKSAANWIJZING, HOUDINGSCORRECTOR
MANUAL DE INSTRUCCIONES, CORRECTOR DE POSTURA
MANUALE D'USO, CORRETTORE DI POSTURA
BRUKSANVISNING, HÅLLNINGSKORREKTUR

PL

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EN

Thank you for purchasing our product. Please carefully read the user manual. Do not use the product before reading and understanding this manual. If you do not understand the warnings, notes, and recommendations, contact a healthcare professional or the seller to avoid product damage or bodily injuries.

WARNING: Check all parts of the product for damage caused during transport. If such damage is found, DO NOT use the product. For more information, contact the seller.

DESCRIPTION AND FEATURES

The posture corrector is made of breathable, lightweight material, ensuring high comfort. It is designed to correct posture disorders and prevent slouching. It helps maintain proper posture by pulling the shoulders back and enforcing correct alignment. Thanks to the Velcro strap system, it can be easily adjusted to the user's individual needs.

CONTRAINDICATIONS

Acute skin inflammations in areas where the posture corrector is in contact with the body (e.g., eczema, ulcers, open wounds, chafing).

Circulation problems – wearing the posture corrector may worsen circulation issues, especially in people with venous or lymphatic disorders.

Osteoporosis – in advanced cases, particularly with already weakened spines, excessive rigid stabilization may lead to fractures.

Respiratory problems – if the brace is too tight or excessively restricts chest movement, it may cause breathing difficulties, which can be dangerous for individuals with asthma, COPD, or other respiratory conditions.

Neurological problems – for patients with nerve damage or sensory disturbances, using the brace may worsen symptoms, e.g., by excessive nerve compression.

APPLICATION

Posture disorders

Back pain syndromes

Healthy spine prevention

SIZE: UNIVERSAL

WARNING:

In case of a “serious incident” related to the product that has directly or indirectly caused, could have caused, or may cause any of the following:

- a) death of the patient, user, or another person, or
- b) temporary or permanent deterioration of health of the patient, user, or another person, or
- c) serious public health risk,

this “serious incident” must be reported to the manufacturer and the relevant authority of the member state where the user or patient resides. In Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices, and Biocidal Products.

WARNING:

If pain, allergic reactions, or other concerning symptoms related to the use of the product occur, consult a healthcare professional.

WARNING: It is prohibited to use the product in any way other than its intended purpose.

WARNING: The manufacturer is not liable for damage caused by neglecting cleaning and maintenance instructions or failure to follow the recommendations in this manual.

TARGET PATIENT GROUP


Individuals dealing with diseases, dysfunctions, or injuries requiring treatment, rehabilitation, or compensation for which this product is intended (see the purpose section in this manual). The product can be purchased independently or based on the recommendation of a doctor, therapist, or another specialist. Whether purchased independently or on recommendation, available sizes, required functions, product variations, indications, and contraindications must be considered, as well as the information provided by the manufacturer.

USAGE

Place the brace on your back and secure it around your waist using the wide Velcro strap.

Next, guide the straps through the back of the brace, pass them through the designated loop, and adjust them by fastening them with Velcro to the front waist area to ensure comfort without causing discomfort.

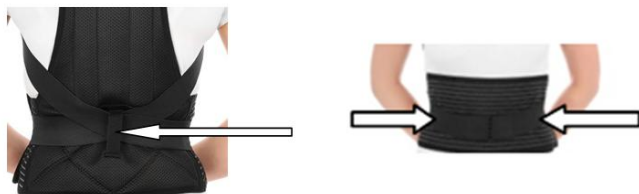You can use additional padding on the shoulders or arms if you feel discomfort or pressure from the straps.

CLEANING AND MAINTENANCE

The product should be hand-washed with a mild detergent at a temperature of up to 40°C. Do not machine wash. Avoid strong detergents, bleach, or abrasive cleaning tools. Dry flat in the shade, avoiding direct sunlight. Do not tumble dry, dry clean, or iron.

STORAGE

Store the product in a cool, dry place, away from moisture and direct sunlight.

DISPOSAL

After use, the medical device can be disposed of as regular household waste

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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