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Model:
Numer Serii:

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR Sp. J., 03-068 Warsaw, ul. Zawiańska 43, hereby grants a 12-month warranty on the rollator, starting from the date of purchase by the Buyer. Products financed by the National Health Fund are covered with the extended warranty period, which covers half of the useful life specified in the Regulation on Guaranteed Services, provided that the useful life is specified in the Regulation in months or years.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, send it immediately to ANTAR (as specified in point 1) or contact the point of sale.
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum permissible weight of the user as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts subject to normal wear and tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED. IN ORDER TO EXTEND WARRANTY PERIOD FOR PRODUCTS FINANCED BY THE NATIONAL HEALTH FUND, ATTACH A COPY OF THE ORDER FOR THE SUPPLY OF MEDICAL DEVICES BEING ORTHOPAEDIC ITEMS AND AIDS.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.

ANTAR Sp. J.
03-068 Warszawa, ul. Zawiańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net
e-mail:antar@antar.net

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Introduction

Thank you for purchasing our product. Before using, please read this user manual carefully. Incorrect use may lead to health detriment, injury or accidents. Therefore, in order to use your wheelchair safely and successfully, please read this manual carefully. Inspect all parts for damage in transit. If any such damage is found, DO NOT use the product. More information from the manufacturer.

Lightweight aluminum IV pole with a mobile base and metal hooks. The tripod has a durable structure that allows easy and quick movement. The stand is characterized by simple height adjustment and the possibility of locking it in the range from 115 to 200 cm. Thanks to the five-star base with wheels, the IV pole is mobile and stable

GENERAL WARNINGS

Before using this product, please read this manual and keep it for future reference. If you are unable to read or understand the warnings, cautions, or instructions, consult your healthcare professional or dealer before using the instrument. Failure to do so may result in injury or damage.

CAUTION: When using and handling the product, and when assembling and adjusting mechanisms, there may be a risk of entrapment and / or pinching of body parts of the user / accompanying person in openings / gaps between the components. These steps should be performed with extreme caution. After completing the adjustment, stabilize the position by tightening the nuts / bolts carefully.

WARNING:

When a product-related "serious incident" occurs, which directly or indirectly led, or could lead to, any of the following events:

- a) death of a patient, user or other person or
 - (b) temporary or permanent deterioration in the health of a patient, user or other person, or
 - c) a serious threat to public health
- the above "serious incident" must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient is resident. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

WARNING:

In the event of pain, allergic reactions or other disturbing, unclear for the user symptoms related to the use of a medical device, consult a healthcare professional.

NOTE: In the event of pain, great discomfort, or an allergic reaction, please contact a healthcare professional as soon as possible.

ATTENTION: The manufacturer is not responsible for damage caused by neglect of maintenance, inadequate servicing, or resulting from failure to comply with the recommendations contained in this manual.

ATTENTION: it is forbidden to use the product other than in accordance with its intended purpose

CAUTION: There may be a risk of tipping over if the product is used improperly. Please follow the boarding / dismounting / moving instructions.

PRODUCT FEATURES:

- light frame, plastic base
- height adjustment in the range of 115 - 200 cm
- mobile base with 5 swivel castors for easy movement
- base diameter 52 cm
- 4 catches made of metal with a max. load 2.5 kg / 1 hook
- frame colors: chrome
- base color: black
- new, medical, certified product
- total weight: 2.58 kg
- maximum payload

This symbol means: "Maximum payload"

DESTINY

- intended for common use in medical centers,
- as hangers for administering all kinds of infusion fluids (so-called drips),
- thanks to the wheels, they can move with the patient

CONTRAINDICATIONS

- physical or mental limitations (such as visual impairment) that prevent you from being safe handling the wheelchair.

STAND ASSEMBLY

There is a hole in the base of the stand to insert the provided screw the adjusting tube is fastened to the base of the stand with a screw

- use the nut to fasten two adjustment tubes
- attach the hook with catches on the top

Height adjustment

- loosen the nut (FIG.1)
 - move the adjusting tube along its axis to the height that suits you
 - check the safety mechanism before next use
- height is properly attached and secured.

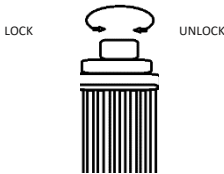
FIG. 1

BEFORE EACH USE, CHECK THAT:

- The locking mechanism is properly installed and secured.
- Wheels are in perfect condition. A worn base plate must be replaced immediately.
- If any damage is found, the product must NOT be used. Failure to do so may result in injury.
- The handle and locking mechanism should not be repeatedly closed and opened unnecessarily. This can handle in place or the locking element or even cause them to break.
- Never try to modify the elbow crutches. They are designed to be safe. For the same reasons, use only the
- The elbow crutches are for walking aid only. Never use them for any other purpose, for example as a lift or

weaken the button that holds the
manufacturer's spare parts.
leverage.

CLEANING AND MAINTENANCE

Failure to care for and maintain the stand can result in serious injury.

- Clean the stand regularly. Wash the handles (handles) and the base (sleeve) with clean water. Wipe all surfaces with a damp, soft cloth and allow them to dry before use. Do not use abrasive products, cleaning agents, chemicals, oily products that could make the balls slippery, or any bleaching products.
- Regularly check the technical condition of all elements of the stand. For safety reasons, the health care professional should carry out an inspection every 6 months. Regularly inspect parts, including hardware, brackets, and plastics for distortion, corrosion, cracks, wear or squeezing.
- Check the wear or breakage of the following parts: base (bushing), handle, locking mechanism. It is important to regularly check the condition of the stand base: depending on the surface on which the stand is used, the stand wears faster or slower.
- For safety and continued use, replace worn or broken parts immediately. Use only manufacturer's replacement parts.
- The only modification you can make to the product is the replacement of worn parts.
- If the product has not been used for an extended period, a healthcare professional should check it again before using it.

STORAGE

Failure to comply with the storage conditions may damage the stand and thus risk of injury.

- Do not store the product in a place directly exposed to extremely high temperatures (eg near a window or near a heater) or low temperatures.
- If any part of the stand is exposed to extreme temperature (above 37 ° C or below 0 ° C), high humidity and / or has been wet before use, check for deformation. Failure to do so may result in damage or injury. Before use, check that all stand components are securely tightened.
- When exposed to extreme temperatures - above 41 ° C - the stand becomes hot and may result in first-degree burns. The ball should cool down before use.
- Store away from fire and flames.
- Comply with the storage conditions for the rack as shown below: - Relative humidity should be between 30% and 70%.
- Protect the stand from dust, corrosion (abrasive elements, sand, sea water, salty air) or an oily environment by storing it in a cover.
- Store in a dry place at a temperature between 5 ° C and 40 ° C.

METHOD OF DISPOSAL OF THE PRODUCT AFTER REMOVING FROM SERVICE

After the device is withdrawn from use, the medical device can be disposed of as normal household waste, except for electrical products - proceed in the manner appropriate for the disposal of electrical and electronic equipment.

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