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TEL: 886-(0)5-2066688 FAX: 886-(0)5-2067788

CE Representative

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28 Newfield Road, Stoubridge, West Midlands DY9 0JR U.K

TEL: +44(0)1562 885097 FAX: +44(0)1562 886838

<http://www.karmamedical.com>

Importer:

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tel. +48 22 518 36 00, +48 22 518 36 15

fax +48 22 518 36 30, +48 22 518 36 31

v1-2 30.09.2022

Manual Wheelchair - User Manual KM-5000

Series with reclining seat

ISO 9001

Lead the way, all the way!

GENERAL WARNINGS

- Keep the truck in good condition, especially check tire pressure, as it affects brake performance.
- Make sure both brakes are applied before sitting on or getting off the wheelchair.
- When the stroller is not being pushed or used, make sure both brakes are applied.
- Use caution when driving on curbs or turning on slopes.
- Avoid steep or uneven surfaces that can cause the cart to tip over.
- Soft surfaces should be avoided, as the front wheels can get stuck in them.
- One should drive calmly regardless of conditions.
- Look ahead to warn of dangers such as furniture, doors, holes, slopes, curbs, front doors, etc.
- Make sure both feet are resting securely on the footrest and clothing is wrapped securely so that it is not pulled in by the wheels.
- Caution should be exercised when driving , downhill and across slopes.
- To control the stroller, the caregiver should keep both hands on the handles to push the stroller.
- Pay attention to the balance of the wheelchair.
- Check the frame and components regularly and use this manual during and after the warranty period.
- If the padding is loose, it should be stretched or replaced with a new one.
- To avoid tipping over, equip the cart with anti-tipping wheels.
- Equip the wheelchair with a handicap safety belt.
- Equip the wheelchair with a calf belt for the disabled.
- Stay away from the fire, as there is a risk of injury to the user and damage to the cart.
- Under no circumstances should the cart be overloaded.
- Do not hang heavy objects from the back of the cart that could upset its stability, especially when climbing an incline.
- Do not allow more than one person to use the cart at the same time.
- Do not use the wheelchair with unconscious patients.
- Do not use another means of transportation to push the cart.
- Do not practice two-wheeled driving, as there is a risk of falling or over and serious injury.
- Do not sit in a wheelchair on a vehicle that is in motion unless it contains restraint hook markings.
- Do not drive or turn too fast.
- Do not drive on a slope with an angle of more than 7 degrees.
- Do not stand on the foot rests when getting on and off wheelchair.
- Do not slide down the steps.
- Do not lift the wheelchair by grasping its moving (removable) parts. Lifting the stroller by its moving (removable) parts may cause injury to the user and/or caregiver or damage to the stroller.
- Do not use the stroller when its rear connector is folded. Folding changes the center of gravity and may cause the stroller to tip over. It should only be used for storing and transporting the stroller.

- Do not hang anything on the rear joint that could cause the backrest to fold down suddenly and cause injury.

MAINTENANCE SERVICES OFFERED BY KARMA

- No tools are required to assemble the cart. The cart should be assembled by an authorized Karma dealer. An authorized Karma dealer also offers maintenance services. Unauthorized repairs, upgrades and additions will void the warranty and may lead to injury. Any rear wheel position adjustment should be performed by an authorized Karma service center. Your Karma service center has the necessary service/maintenance instructions for this model of wheelchair.
- In case of defects in the parts of the wheelchair or their loss, contact the dealer for repair and/or replacement. Use only parts authorized by Karma. Use of parts that are not authorized by Karma will void the warranty.
- Customers should order parts, repairs and periodic inspections from their dealer.
- To avoid voiding the warranty, customers should perform all service work at their dealer. Always contact an authorized service center. Unauthorized repairs will void the warranty. Any procedures not described in this manual must be carried out by a qualified professional. He, too, has a service manual with a complete list of maintenance steps required to keep your particular Karma stroller model in good condition. For information regarding Karma's warranty and payment for maintenance and parts, see the "Warranty" section.
- At moment, Karma does not have an official list of service points. For service matters, please contact your dealer. If this is not possible, you should contact another Karma-authorized dealer for service matters. If this is not possible, contact directly the company that will arrange for service repair. Note that there will be a charge for pick-up and delivery of the cart.

CE SIGN

This cart complies with the requirements of Directive 93/42/EEC on medical devices.

NOTES:

Before purchasing and using a Karma wheelchair, a qualified expert should be consulted to ensure that the product chosen is suitable and that training on its use will be provided. It may be necessary to equip the cart with additional parts or change its settings. Consult your service center on this matter.

TABLE OF CONTENTS

1. INTRODUCTION.....	6
2. SECURITY	6
3. USE AS INTENDED AND PARTS OF	10
4. USE.....	13
5. INSPECTION AND MAINTENANCE.....	25
6. OPTIONS AND ACCESSORIES.....	28
7. PROBLEM SOLVING	30
8. SPECIFICATIONS.....	31
9. CLEANING AND RECYCLING	31
10. WARRANTY	33

1. INTRODUCTION

Read these instructions carefully before using the wheelchair. Incorrect use of the stroller may result in personal injury. Therefore, in order to use it in a carefree and pleasant manner, please read this manual.

- This manual contains operating instructions for the wheelchair, tips for and folding, and for dealing with possible accidents.
- This manual applies to the Karma model KM-5000 wheelchair with reclining seat.
- The symbols used in this manual are explained below. Pay special attention to the parts of the manual containing the described symbols.

	WARNING:	Incorrect use can cause serious injury or death.
	NOTES:	Incorrect use may cause serious injury and/or damage to the wheelchair.
	SUGGESTION:	Follow these tips to keep your cart in good condition

- This manual contains a table of repair and maintenance activities as well as details of the guar-
- If another person is using the wheelchair, make sure they have access to this manual.
- Since stroller models are subject to change, some of the illustrations and pictures appearing in manual may not correspond to the design of the stroller purchased. The manufacturer reserves the right to make changes to models without prior notice. This manual will be updated regularly, the latest version is available on Karma's website.

2. SECURITY

2.1 Door It is safest to use the wheelchair on level surfaces that do not contain obstacles. If possible, use a ramp to enter the premises. Keeping your hands on the push handles allow you to maintain maximum balance. Do not overcome obstacles by force, by means of propulsion, and under no circumstances use the sides of the entrance door as an aid to pull yourself through the entrance.

2.2 Leaning out

Do not lean out of the wheelchair, as this may cause you to lose your balance and fall (see Figure 2.1).

2.3 Control at inclines

- Never drive up a hill, due to risk of tipping.
- When pushing the wheelchair, go straight up the slope. If you need to descend, drive the cart sideways until the cart user is on a flat surface (see Figures 2.2. and 2.3)
- Avoid stopping on sloping terrain, as you may risk losing control of the cart.
- Do not use the parking brakes to slow or stop the forklift; this may provide the forklift to seize or tilt to one side, resulting in loss of control.
- Whenever possible, in sloping areas, steer the wheelchair with the help of a caregiver.

2.4 Curbs

- Whenever possible, use an attendant/assistant (if using an attendant, be sure the wheelchair is equipped with push handles).
- When driving off curbs, avoid hard impacts. A strong impact can damage cart.
- Under no circumstances should you drive on or off curbs higher than standard curbs (30-50 mm).
- Always use caution when stepping onto the curb and do not back up or down the curb without assistance from a caregiver.

2.5 Stairs

- Always use two or more caregivers. The wheelchair should be placed in balance. One of the caregivers (standing at the back) supports the wheelchair at the first step, firmly holding the its handrails. The second caregiver, holding the stationary part of the front frame firmly, lifts the wheelchair over the stairs and holds it while the first caregiver puts one foot on the next step and repeats his operation (see figure below).
- Never go down the stairs backwards.
- Never use an escalator; use an elevator instead.

2.6 Streets

- Always use extra caution when crossing the street. Reflective tape affixed to the stroller and/or the user's clothing will help ensure that the stroller is more visible to drivers.
- Avoid potholes and uneven areas, which will help cause the front wheels to jam (see Fig. 2.4)

2.7 General comments


- Always make sure that BOTH brakes are applied before sitting on and getting off the wheelchair.
- Under no circumstances should you stand on the foot rests when getting on and off wheelchair.
- Under no circumstances should more than one person be allowed to sit in the stroller. Karma strollers are designed for a single user.
- Make sure that clothing does not get caught in the rear wheels.

2.8 Cladding

- Use warm water with mild soap to clean the upholstery.
- You can also use upholstery cleaners/foams, but avoid other cleaning solvents.

2.9 Casters

anti-tip wheels You should purchase anti-tip wheels and use them every time to avoid injury. The distance between the rear wheel and the anti-tip wheel should not be less than 5 cm. The anti-tip wheel should not be fixed higher than 5 cm above the ground.

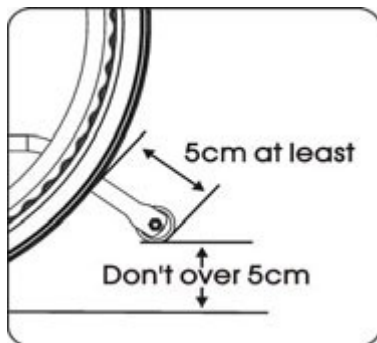2.10 Maximum permissible user weight

The KM-5000 series is suitable for users weighing no more than 115 kg. Please be that loads exceeding the weight limit may damage the wheelchair and cause malfunctions, creating hazards. The warranty does not cover damage caused by improper use of the wheelchair.

2.11 General usage errors should be avoided:

- Keep in mind that any vehicle used incorrectly can cause injury. The following manual contains numerous warnings. To increase safety, proceed in accordance with the guidelines. Be aware that unreasonable use endangers the safety of the user and his surroundings.
- In open spaces, always observe the Traffic Law and the guidelines in this manual. Used carefully and judiciously, the wheelchair will serve its user well. k for many years to come.
- Do not allow children to stand on or play with the stroller (see Figure 2.5).

Figure 2.1

Figure 2.2

Figure 2.3

Figure 2.4

Figure 2.5

2.12 Marking

Before using the wheelchair, carefully read all its markings. Do not remove them, but keep them for future reference.

Before

- The rear spoke wheels are about 70 cm in diameter and have tires made of microcellular rubber.
- For use of the wheelchair by full or partial amputees (above or one knee or both knees) or for other circumstances affecting the users' natural center of gravity, it is suggested to contact a therapist or retailer as well as to attach anti-tip wheels. The wheelchair is designed so that the rear wheels can be reversed for amputee users.

3. USE AS INTENDED AND PARTS

3.1 Use as intended

- The KM-5000 series includes self-propelled and caregiver-driven wheelchairs. This means that they can be driven by the caregiver or by the user (self-propelled).
- The KM-5000 series is suitable for users weighing up to 115 kg. The tilt of the backrest can be adjusted from 90°-163°.
- The KM-5000 series is an adjustable cart. It is suitable for users who require side-to-side transfer. The moving parts are easy to operate.
- It has solid tires. It also has a frame made of extremely durable aluminum alloy and hard accessories that are resistant to extreme temperature fluctuations (neither the user nor the caregiver touch the metal parts during normal use).
- The front wheels have polyurethane tires for greater durability; pneumatic wheels are available upon request.

3.2 Part description

1. Pillow
2. Headrest
3. Backrest padding
4. Push handle
5. Folding lever
6. Armrest
7. Parking brake
8. Rear wheel

9. Anti-tip wheel
10. Frame/background
11. Front wheel
12. Seat padding
13. Footrest lifting lever
14. Calf support
15. Foot rest

- Rear wheels with a radius of 14 cm available.

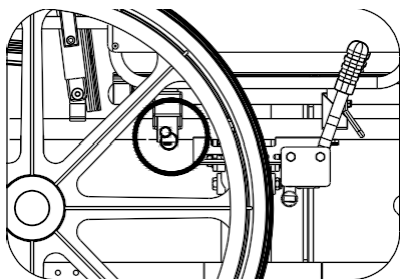
Figure 3.2 Self-closing gas spring

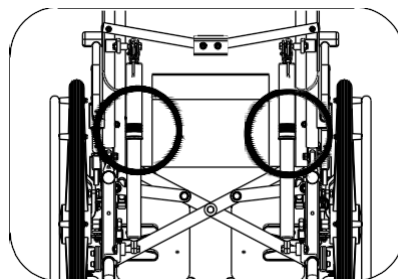
Figure 3.3 Swiveling circular latch

4. USE

4.1 Preparation for use

4.1.1 Unfolding

- Supporting the armrests, pull them outward (see Fig. 4.1)
- Pointing your fingers inward and keeping the center of your hand on the frame, press down until the seat fully lowers. Make sure that the frame is rigidly placed in the part that holds it. shed (see Fig. 4.2)
- Fix the footrest and check that it is locked (see Fig. 4.3)
- Fold out the foot rests (see Figure 4.4)
- Press on the transverse supporting part (see Fig. 4.5)
- The cart is ready for use! (see fig. 4.6) Alternatively, you can stand next to the cart and take step one by pulling closer to you Armrest and push back further.

Figure 4.1

Figure 4.2

Figure 4.3

Figure 4.4

Figure 4.5

Figure 4.6

NOTES:

When performing step 2, under no circumstances should you hold your fingers on the sides of the seat or under it, as they may be caught between the seat and the frame (see Fig. 4.7)

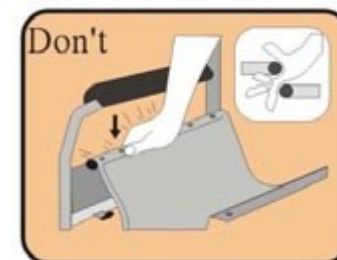
Figure 4.7

4.1.2 Headrest assembly

- Put the headrest tubes on the tail pipes (see Fig. 4.8)
- Press the pipe locking button (see Fig. 4.9)

Figure 4.8

Figure 4.9

4.1.3 Filing

- Fold the foot rests (see Fig. 4.10)
- Roll up and remove the footrest (see Fig. 4.11)
- The transverse supporting part should be loosened (see Fig. 4.12)
- Remove the headrest (see Fig. 4.13)
- Pull up the seat padding by holding it by the center or by the front and back edges (see Fig. 4.14)
- Folding is complete (see Fig. 4.15)

Figure 4.10

Figure 4.11

Figure 4.12

Figure 4.13

Figure 4.14

Figure 4.15

NOTES:

Be careful to prevent your fingers from getting caught when folding the stroller.

4.1.4 Adjustment

4.1.4.1 Armrests

- Standard armrest: removable armrest
- The armrest can be removed by pulling out the self-locking latch and releasing the lever under the armrest. The armrest can then be removed by lifting it up (see fig. 4.16 i 4.17)

Figure 4.16

Figure 4.17

NOTES:

Before using the armrest, make sure that the self-closing circular latch is positioned in the hole in the armrest (see Fig. 3.2)

4.1.4.2 Footrest lift

- You can raise the footrest in the vertical direction and set it to any desired height (see Fig. 4.18)
- If need to lower the footrest, press the release lever (see Fig. 4.19)
- To adjust the height of the footrests, use Karma's tools to loosen the bolt on the side of the footrest.
- Adjust them to the required height and lock them again (see Fig. 4.20)

SUGGESTION:

Make sure that the height of the footrest is properly adjusted and corresponds the length of the user's calf, so as to reduce pressure on the body.
The footrest must be at a height higher than 5 cm (see Fig. 4.21)

Figure 4.18

Figure 4.19

Figure 4.20

Figure 4.21

NOTES:

When adjusting the height of the footrest, make sure that the user's thighs are safely the seat and that the inside of the knee does not hit the front edge of the seat (see Figures 4-22 and 4-23)

Figure 4.22

Figure 4.23

Figure 4.24

4.1.4.3 Installation of an anti-tip wheel

Refer to the figure below and make sure that the button is locked into the tube. Also test the operation of the wheels (see Figure 4.24)

4.1.3.4 Rear wheel

The cart can be equipped with wheels on a quick-release axle. In this case, follow the instructions below.

4.1.3.4.1 Removing the rear wheel

- Balance the weight of the stroller with one hand while holding the push handle.
- Press the wheel hub cover and unlock the quick release axle mechanism.
- Remove the wheel (see Fig. 4.25)

4.1.3.4.2 Rear wheel overlay

- First, make sure the brake is in the OFF position. Raise the wheel with the axle and press the button in the center of the axle as above. This will retract the ball bearings, allowing you to push the axle into the wheel bracket located at the bottom of the rear frame of the wheelchair. When removing the thumb from the axle button, the ball bearings should visibly protrude beyond the axle at the rear of the wheel bracket. This will allow the wheel to be secured in the frame (see Fig. 4.26)

Figure 4.25

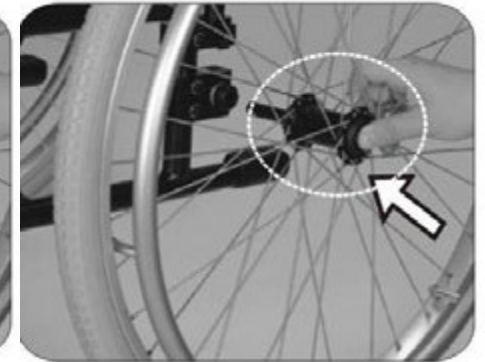
Figure 4.26

4.1.3.5 Important recommendations and prohibitions

- You should: keep axles clean and oiled with light grease or maintenance spray
- You should: maintain adequate tire pressure for pneumatic tires
- Do not: attempt to remove the wheels when the wheelchair is in use, even if the over-turn wheels are used.

4.1.3.6 Seat belt

How to install and remove a seat belt

- Loosen the fourth and fifth screws on either side of the seat padding and insert the seat belt between the fourth and fifth screw holes (see Fig. 4.27)
- For comfortable use, fasten the belt at a 45° angle (see Fig. 4.28)
- Put the seat back on and tighten the screws.

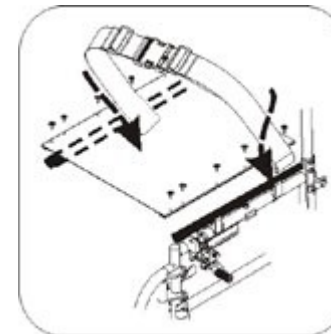
Figure 4.27

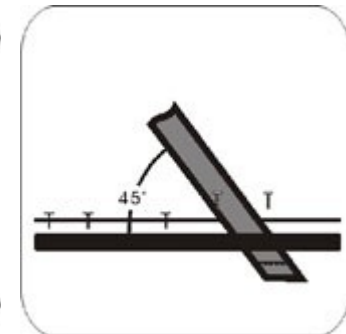
Figure 4.28

Seat belt length adjustment

- Specify the required length of the belt (Fig. 4.29)
- To lengthen or shorten the belt, slide the double part of the belt through the buckle (Fig. 4.30)
- Pull the double part of the belt firmly (Fig. 4.31)
- Pull the buckle up to the plastic clasp so that the belt cannot slip out of its position.

NOTES:

- Make sure that the axle is fully seated in the wheel bracket and that the ball bearings are fully open and have properly secured the wheel in the frame.
- The rear wheel can be fixed in the frame in different positions according to the height of the seat and for the overall stability of the wheelchair. Changing the position of the wheel should be done
A mechanic or therapist with the appropriate qualifications.

Make sure that there is no space between the belt and the buckle if the wheelchair user is carried by different people (see Fig. 4.32)

SUGGESTION:

It is recommended that the belt be secured quite tightly around the hips. It should be tightened until there is space between the wearer and the belt to allow a finger to be inserted. Make sure the belt is properly fastened recommended. If this is not possible, contact your dealer.

Figure 4.29

Figure 4.30

Figure 4.31

Figure 4.32

How to use the belt:

- Be sure to fasten the seat belt while sitting in the stroller.
- Adjust the belt and fit it the body ensure comfort and safety.
- When not in use, it should be fastened so that it does not screw into the rear wheels when moving.

Maintenance

The seat belt is made of 100% polyester fiber and its fastener is made of plastic. It requires no special maintenance, but should be checked monthly to make sure its edges are not frayed and weekly to make sure it fits the wearer. If there are any signs of damage, it should be replaced immediately.

4.1.4 Transportation and storage

sure the cart is folded when transporting it

SUGGESTION:

Store the cart in a place out of sunlight, and not exposed to rain or dew. Consult your Karma dealer for details.

WARNING:

Do not place objects on the cart when it is stored. Holding objects on a folded or unfolded cart will cause damage to the cart that is not covered by the warranty (see illustration at right)

4.2 Putting a person in and out of a wheelchair

4.2.1 Putting a person in a wheelchair

- Both handbrakes should be applied.
- Turn the foot rests upside down and twist them to the side.
- Keep the user low over the seat so that they can support themselves on the armrests.
- Turn the foot rests back down so that the user can place their feet on them. Make sure the heels are resting on the backrests to avoid contact with the front wheels.
- Both handbrakes should be released.

NOTES:

Always make sure that the front wheels are aligned with the rear wheels and pointing forward, and that the back of the seat is belayed by a wall or held firmly by the caregiver.

4.2.2 Removing a person from a wheelchair

- Make sure the cart is close to the user's destination.
- Position the wheelchair backwards so that the front wheels face forward and align with the rear wheels. Ensure that, if possible, the back of the seat rests against the wall a If this is not possible, it must be held firmly by the caregiver.
- Both handbrakes should be applied.
- Turn the foot rests upside down and twist them to the side.
- You should lower yourself onto the seat and place one foot firmly on the ground while maintaining the proper angle of knee flexion, and place the other foot at a small distance and slightly behind.
- Place your hands on the front of the armrests and, shifting your body weight forward, pull your head and shoulders forward over your knees and lift yourself up.

Some people require getting in and out sideways.

- Tilt the armrests back and make sure the wheelchair stands close to its destination.
- Support the body and slowly move sideways, using the board to move. Be careful not to hit the brakes or any other part of the wheelchair.

NOTES:

- Never stand on the foot rests when getting on or off the wheelchair (see Figure 4.33)
- Use handbrake when getting on and off the wheelchair. The hand brake is applied using the lever on the side of the wheelchair.
- Make sure that neither tire treads nor brake pads have worn out. If pneumatic tires are used, to avoid premature wear, you should maintain the recommended pressure.
- Make sure that the user's fingers not caught by the brake.

To prevent the wheelchair from rolling away when getting on and off the wheelchair, make the brakes are applied. To ensure easy ascent and descent, make sure the foot rests are turned upward (see Figure 4.34)

4.3 General use

4.3.1 Caregiver push vs. autonomous wheelchair drive

- Make sure that the brakes are locked whenever the wheelchair is not in use (see Fig. 4.35)
- Always make sure that both feet are resting securely on the foot rests, and that clothing is safely rolled up and in no danger of being pulled into the wheels (see Figure 4.36)
- To control the stroller, always hold the push handles with both hands (see Fig. 4.37)
- Steep or uneven surfaces should be avoided, as there is a risk of the cart tipping over (see Fig. 4.38)
- Soft surfaces should be avoided, as the front wheels may get stuck in it (see Fig. 4.39)
- Always be cautious and do not move or turn too fast, trying to maintain a steady and even pace of movement.
- Before crossing the street, stop completely and look both ways.
- Drive through ruts at an angle, which will prevent the front wheels from locking.
- Do not cross railroad tracks at full speed.
- Keep your hands and arms inside the wheelchair when passing through the front door.

THEN maintain contact with the person in the . When a caregiver stops to look at exhibits or talk, he or she should remember to turn the wheelchair so that the person in the wheelchair does not have to turn around.

WARNING:

Under no circumstances should the brakes be applied too abruptly, as this may cause the user to be thrown from the wheelchair.

Figure 4.33

Figure 4.34

Figure 4.35

Figure 4.33

Figure 4.36

Figure 4.34

Figure 4.37

Figure 4.35

Figure 4.38

Figure 4.39

Figure 4.40

4.3.2 Tilting the wheelchair by the caregiver

- Make sure that both of the user's feet are on the foot rests and do not slip off, and that the hands have a firm grip and do not slip off the handles.
- Both brakes should be released.
- Both push handles should be grasped firmly.
- It is necessary to put the stronger foot on the appropriate tipping tube at the back of the wheelchair.
- It is necessary to pull the cart backward with both hands and press downward and forward with the foot. This will put the cart on the rear wheels (see Fig. 4.40)

NOTES:

- When the stroller is tilted at an angle of about 25°, it will keep its balance on the rear wheels and be easy to hold. If it is tilted at a higher angle, It will be difficult to hold it, and the person in the wheelchair will feel insecure.
- You should lower the wheelchair with your back straight while keeping your foot on the anti-roll wheel and lowering the wheelchair slowly and smoothly to the ground.

4.3.3 Tilting the cart by the user

- Make sure your feet rest on foot rests and do not slip off.
- Both brakes should be released.
- Grasp both push rails and push resolutely forward directly in front of the obstacle and continue the pace of pushing (you should be ready to push continuously) until moment-Here, when the obstacle will be overcome.

NOTES:

- The above method should only be used by people experienced in using a wheelchair. Strong arms are necessary. It is always better to use caregiver assistance.
- Leaning backwards causes the front wheels to lift BUT is much safer than leaning forward and maintaining balance.

4.3.4 Descending from

curbs Note that it is always better, if possible, to avoid curbs by using ramps or , and the following operation is best carried out with the help of a caregiver.

4.3.4.1 Sliding off curbs ahead with the help of a caregiver

NOTES:

- This method is more difficult and can cause jerkiness, but it allows a better view of the traffic situation.
- Care should be taken to ensure that the anti-tip wheel does not lock on the curb causing the user to fall.

UNDER NO CIRCUMSTANCES SHOULD THE WHEELCHAIR BE ALLOWED TO TIP FORWARD WHEN STEPPING OFF THE CURB.

- Position the wheelchair so that the front wheels are close to the edge of the curb and make sure that the road is clear.
- Move the wheelchair slowly and carefully forward, while tilting it slightly backward. Keep your weight ready when the rear wheels cross the edge of the curb. To avoid tipping over, make sure that both rear wheels touch the ground at the same time.
- Gently lower the front wheels. JERKING SHOULD BE AVOIDED.

4.3.4.2 Sliding off curbs in reverse with the help of a caregiver

NOTES:

This method is easier and gentler BUT requires the caregiver to step backwards onto the road.

- Turn the wheelchair with its back to the edge of the curb and ensure that the rear wheels are from the edge.
- While grasping the push handles, gently lower the rear wheels to the road. Be ready to control the weight when the stroller lowers, and ensure that the rear wheels touch the ground in the at the same time.
- Tilt the stroller slightly backward and pull it onto the road until the front wheels and footrest could be lowered freely from the curb to the ground.

4.3.5 Driving onto curbs Note that this operation is best carried out with the help of a caregiver, and the best method of driving onto curbs is forward driving, as this allows you to maintain the direction of movement and is less physically demanding for the caregiver.

4.3.5.1 Driving onto curbs ahead with the help of a caregiver

- Tilt the stroller at the point where the footrest almost touches the curb. The stroller should be tilted to a height that allows the front wheels to drive freely onto the sidewalk.
- Push the stroller forward until the rear wheels touch the curb. Then continue the movement and, while pushing the stroller forward, lift it by the push handles.

4.3.5.2 Riding up curbs in reverse with the help of a caregiver

It should be noted that this method is not recommended for safety reasons.

4.3.6 Riding up and down stairs/stairs

- Stairs should be avoided especially if assistance is not available.
- If possible, use ramps or elevators.
- Two caregivers are needed, each able to carry at least half the combined weight of the person and the wheelchair.

- Each caregiver should bend his knees and grasp the main frame: caregivers in the front should keep their hands directly under the front of the seat, and those in the back should hold the frame directly, nio under the push handles.

NOTES:

Do not lift the wheelchair by grasping the moving (removable) part. Lifting the stroller in this manner may cause injury to the user and/or caregiver or damage to the stroller.

LIFT THE WHEELCHAIR WITH AN UPRIGHT BACK

To lift the wheelchair, straighten your legs. Carry it or carry it up the stairs carefully paying attention to foot placement - do not cross your feet to avoid tripping.

4.4 Transport

The KM-5000 series has passed crash tests in accordance with ISO 7176-19 (Report No. MBK 02/1040) and is suitable for use in vehicles that have been adapted for this purpose. The truck must be securely as described in this manual. The KM-5000 series has been tested using a four-point attachment system. Refer to the manufacturer's manual for more information on securing the cart.

Karma always recommends that the user be carried to the seat of the car with a belt or harness fastened while driving.

SUGGESTION:

possible, move the user to a car seat.

4.5.1 Before using

NOTES:

The KM-5000 series has been tested using a four-point attachment system with the back of the wheelchair in an upright position. Leave the headrest in place during transport.

- Remove all easily removable accessories of the stroller and put them in a safe place where they will not be in motion, so that they do not pose a danger.
- The four-point attachment system should be clipped tightly to the bottom of the wheelchair frame marked with the yellow sticker shown in the picture.

- To ensure user safety, the user's lap belts must also be secured to the transport vehicle. Follow the instructions for use of the lashing system-
nia.
- To ensure the safety of the wheelchair user, the must recheck all fastening points for proper connectivity and sufficient tightness.
The caregiver/transporter must be fully competent in operating the fastening system and aware of its weight limitations.
- Further information available at <http://www.unwin-safety.com/>

4.5.2 Attachment method


The KM-5000 series has been tested using four-point attachment systems with the rear of the wheelchair in an upright position. Details on this operation are available in the system manufacturer's manual.

Putting on the front mount

- You need to connect the straps around the frame of the wheelchair and the front wheel support, as in the picture below. The attachment point is marked with a yellow sticker.
- Remove the brakes and rotate the cart to tighten the front straps. Then apply the brakes (during transport, the brakes should be applied)


Putting on the rear mount

- Attach the back straps to the rail about 300 mm behind the back of the wheelchair.
- On the back of the cart, hooks should be attached to the bracket as shown in the picture. Be sure to stretch the attachment so that the cart is secured.

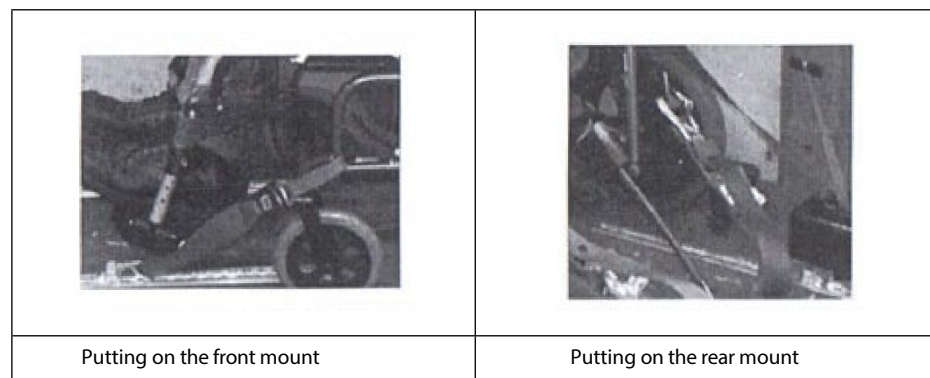Putting on the seat belt

- A three-point seat belt should be fastened around the user.
- The wheelchair in which the person is seated should be facing forward with the parking brake applied.
- The hip mounts should touch the front of the body with the entire surface until the thigh and hip meet, as shown in the image below (correct position of the seat belt)
- Attachments should not be separated from the body by parts of the wheelchair, such as armrests or wheels, as shown in the image below (incorrect seat belt position)

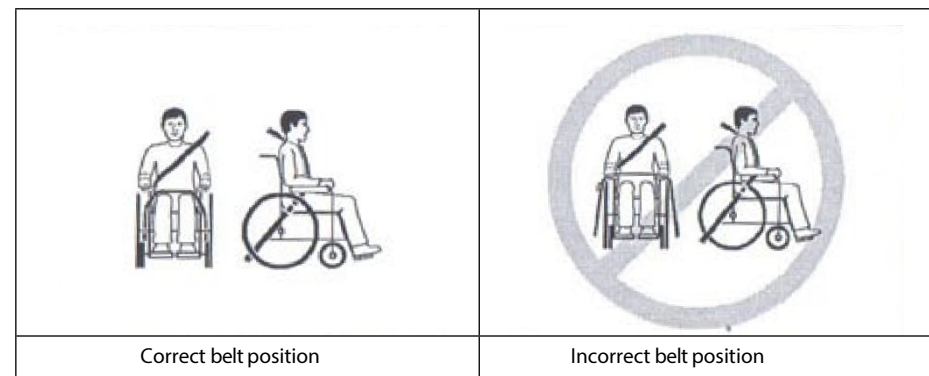5. INSPECTION AND MAINTENANCE

5.1 General

Before traveling, inspect the parts described below. If any abnormality is found, visit your Karma dealer for additional inspection before use.

- Check the brakes for wear, damage or misalignment. Replacement or adjustment of brakes should only be performed by a qualified mechanic.
- Check the hand grips for damage and the attachments on the handles.
- Check the anti-tip wheels damage and secure attachment.
- Check that the anti-tip wheels and frame are not damaged and are securely installed.
- If the cart has been used for a short period of time, its user should check that the bolts and bolts are tightened. All parts should be checked for damage and wear, and if co- If necessary, replace them with the help of a qualified specialist. Check all parts for fit.
- Check the brakes for proper operation.
- Wheels and tires should be checked periodically for cracks and wear and replaced by a qualified professional.
- Qualified specialists are available from an authorized representative.

NOTES:

If any abnormalities are found, visit your Karma cart dealer for inspection and maintenance.

5.2 Documentation regular activities maintenance To make sure your cart is in good condition, visit your Karma dealer regularly for maintenance (and keep records of it) every six months after purchase. Service/maintenance may be subject to a fee.

SUGGESTION:

Even if the cart is not used for a long time, you still need to perform annual maintenance.

Part of	Weekly inspection	Monthly check
GENERAL		
• Does the wheelchair go straight (no erroneous tipping or turning to one side)		V
FRAME AND CROSS MEMBER		
• Checking the frame structure for small cracks/crevices		V
• Checking the bending of the frame or cross member		V
BRAKES		
• Do they not go into the tires while moving		V
• That the pressure pins are not worn or loose		V
• Are they easy to pull on	V	
• When they are pulled, do they prevent the wheelchair from moving	V	
SEAT AND BACKREST		
• Check for rippling and collapse	V	
• Check for a loose or broken base, as a frame, etc.		V
FRONT WHEELS		
• Controlling the smoothness of spinning		V
• Checking for worn bearings or mounts and worn or damaged tires		V
REAR WHEELS		
• Does the quick release axle lock properly; oil if necessary	V	
• Is there no side-to-side twisting when lifting and turning or deviation		V
• Check for cracks, bends or broken spokes	V	
ARMBOARDS		
• Control of the effective locking of the armrest release button		V
• Making sure the armrest detaches and folds properly backward		V
FOOTREST		
• Checking the proper tightening of the screws used to adjust the height of the foot rests		V
• Checking the correct locking of the footrest release lever		V
TIRES		
• Inspection of air loss points, cracks and wear	V	
• When using pneumatic tires, check for proper inflation	V	
• Check if the tread depth is more than 1 mm	V	
• Tire pressure control, pressure values are located on the side of the tire	V	V

NOTES:

Never disassemble or alter the frame structure of the wheelchair. Doing so risks injury or damage. Doing so will void the warranty.

6. OPTIONS AND ACCESSORIES

	
▪ Magnesium alloy rear wheel (MG-40) ▪ Silver spokes+ bright yellow opo-on. ▪ Fits all Karma hand trucks	▪ Drip stand ▪ Height can be adjusted from 130 cm to 200 cm ▪ Fits all Karma hand trucks
	
▪ Oxygen tank ▪ Standard size fits behind cart ▪ Fits all Karma hand trucks	▪ Backrest bag ▪ Standard size ▪ Fits all Karma hand trucks
	

	
▪ Collapsible beverage container ▪ It folds and adjusts to the size of the bu- telcos (6-10cm) ▪ Fits all Karma hand trucks	▪ Castor anti-turn wheel Fits all Karma hand trucks
	
▪ Multifunctional seat belt Fits all hand trucks Karma companies	

NOTES:

Putting on accessories can disturb the overall stability and center of gravity of the wheelchair. Therefore, accessories should be put on by a qualified mechanic and therapist.

7. SOLVING PROBLEMS

Problems related to the use of the cart

Before returning the cart to the dealer, you can try the following solutions to problems.

Problem	Checkpoints and corrective actions
1. cart turns right/left	• Checking tires for proper and even pressure • Checking the angle of the front wheel support • Check that the rear wheels are at an equal distance from the frame seats
2. slow movement or twisting	• Checking tires for proper and even pressure • Check for loose nuts and bolts
3. wobbling/wobbling front wheels	• Checking tires for proper and even pressure • Check for loose nuts and bolts • Checking the angle of the front wheel support
4. creaking and knocking	• Check for loose nuts and bolts
5. loose seat	• Check for loose nuts and bolts

If you have been in solving the problem, contact your cart dealer for assistance.

8. SPECIFICATIONS

Model	KM-5000 F24"
Total length (mm/inch)	1265 mm (49.9")
Total width (mm/inch)	650 mm (25.6")
Overall height (mm/inch)	1290mm (50.8") (930mm/36.6" without headrest)
Total weight (kg/pounds)	22.1 kg (48.7 lb)
Seat width (mm/inch)	41, 43, 46 cm
Seat depth (mm/inch)	554 mm (17.5")
Seat height (mm/inch)	505 mm (19.9")
Headrest height (mm/inch)	455 mm (17.5")
Armrest height (mm/inch)	230 mm (9.1")
Calf length (mm/inch)	410 mm (16.1")~485 mm (19.")
Total length after folding (mm/inch)	910 mm (35.8")
Total width when folded (mm/inch)	390 mm (15.4")
Total height when folded (mm/inch)	930 mm (36.6")
Maximum allowable user weight (kg/lb)	115kg (255lb)

Static stability (According to ISO 7176-1:1999)

Front	Max. 25	Min. 21 degrees
Rear	Max. 24	Min. 22 degrees
Side	Max. 20	Min. 21 degrees

- The degree tolerance is +0, -1 degrees
- Test conditions: the weight of the model used for the test - 115 kg.

NOTES:

Do not drive across a slope (in any direction) with a slope angle than 7 degrees.

9. CLEANING AND RECYCLING

9.1 CLEANING

In order to extend the life of the stroller, it is suggested to clean it periodically (if not daily), especially after exposing it to rain or snow. It is recommended to avoid using the stroller during rain or snow. If the user is surprised by rain while traveling, he should try to find shelter. If this is not possible, continue the journey to the destination point.

Although the frame of the stroller is aluminum, some parts of the frame such as bolts, nuts and wheel spokes will rust after exposure to water. After exposing the stroller to water, dry it very thoroughly. Be sure to store it in a dry place. Do not clean it with abrasive materials.

9.1.1 Frame

Wipe with a damp cloth.

9.1.2 Cladding

- Use warm water with mild soap to clean the upholstery.
- You can also use an upholstery cleaner/foam, but avoid cleaning .

9.2 RECYCLING

The recycling plan should comply with Environmental Protection Agency regulations. The items in the following list are for informational purposes only, as they may differ from local laws and regulations.

Material	Name of Part	Method
Aluminum alloy	Frame	Recyclable
Plastic	Modular elements	Recyclable
Paper	User Manual	Recyclable

- There are many organizations in many countries, knowledge of which is available on the Internet, that accept used goods for charity. You should take advantage of the infostrada's capabilities and help those in less favorable location, protecting our planet for future generations. If you don't have access to the Internet, you can use a local phone line operator or get information from a government center. Earth Day is every day.
- Check with your local recycling center for information on where and how to recycle your cart, as details may vary by region.

10. GUARANTEE

10.1 Number

serial To ensure service and warranty after the sale of the cart, fill in the field by writing the serial number located on the cart frame and keep the entry.

Model	KM-5000 series
Serial number	

10.2 Warranty policy

10.2.1 Content and duration of the guarantee

- The repair is free of charge if the problems are due to manufacturing defects in the parts and occur within warranty period (this type of repair will be referred to as "warranty repair").
- Warranty repairs can only be carried out by authorized Karma dealers.
- Karma is not responsible for the cost of transporting the cart repair or replacement.
- The parts listed are the property of Karma.

10.2.2 Problems not covered by the warranty

- Irregularities that do not affect the operation and safety of the product.
- Damage resulting from normal wear and tear.
- Problems arising from the lack of basic maintenance and self-inspection of the truck.
- Damage resulting from failure to follow usage guidelines.
- Installation or use of parts not manufactured by Karma Company.
- Damage resulting from use under unsuitable conditions (see the guidelines of this manual).
- Any additions, modifications or changes not authorized by Karma.
- Damage from external agents such as chemicals, oil, debris, etc.
- Natural disasters.

10.2.3 Karma does not raise responsibility for the following costs

- Replacement or replenishment of consumables, grease, fluid.
- Inspection, adjustment, cleaning, etc.
- Costs of self-inspection recommended by Karma.

10.2.4 Responsibilities

10.2.4.1 User responsibilities

- Fully comply with the usage guidelines in this manual.
- Perform self-inspections and basic maintenance.
- Keep records of self-inspections and basic maintenance.
- Do not allow any changes, additions or modifications or use of parts not manufactured by Karma.
- Carefully read and understand the instructions.

10.2.4.2 Responsibilities of the seller

- Follow the guidelines of this manual and provide the user with information on product and basic maintenance.
- Verify that the warranty conditions have been met.
- Provide relevant warranty information and recommendations for improvements.

10.2.5 Disclaimer

Parts that easily, such as tires, front rollers, seat cushions, backrest cushions, etc., are not covered by the warranty.

10.2.6 Duration of warranty and parts covered

your local dealer.

10.2.7 Transfer of warranty rights

If you purchase a used Karma stroller that is still under warranty, you must register the stroller with the dealer where it was purchased, based on the presentation of the stroller and this manual.

:	
Gender:	
Date of birth:	
Address:	
Model:	
Serial number:	
Date of Purchase	

Vendor details:

Store Name:	
Phone and address:	

If you have suggestions on how we can improve our products, we encourage you contact us and share your feedback on the wheelchair you purchased. Our contact information is available at www.karmamobility.co.uk. Thank you and we wish you pleasant use.

KARMA MEDICAL PRODUCTS CO.,LTD.

Number 2363 section. 2 Da-Shiue Rd., Min-Hsiung Shiang, Chia-Yi, Taiwan

PHONE: 886-(0)5-206 6688 FAX: 886-(0)5-206 7788

CE Representative

28 Newfield Road, Stoubridge, West Midlands DY9 OJR U.K

TEL: +44(0)1562 885097 FAX: +44(0)1562 886838

<http://www.karmamedical.com>

CONTRAINDICATIONS

physical or mental limitations (e.g., visual impairment) that prevent safe handling.

CAUTION: If the product is not used properly, there may be a risk of over. Please follow the instructions for getting on/sitting down/moving.

CAUTION: During the use and operation of the product and when assembling and adjusting the mechanisms, there may be a danger of trapping and/or squeezing parts of the user's/companion's body in the openings/gaps between components. Perform these operations with particular care. Stabilize the position by carefully tightening the nuts/bolts after the adjustment is completed.

CAUTION: it is forbidden to use the product other than for its intended purpose.

NOTE: Before use, check that all components have been assembled correctly.

NOTE: The manufacturer is not responsible for damage caused by negligent maintenance, inadequate servicing or resulting from failure to follow recommendations in this manual.

NOTES:

In the event of a product-related "serious " that directly or indirectly led to, may have led to, or may lead to any of the following:

a) death of a patient, user or other person or

b) temporary or permanent deterioration in the health of a patient, user or other person or

c) a serious threat to public health

the above "serious incident" should be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the case of Poland, the competent authority is the Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

NOTES:

If you experience pain, allergic reactions, or other disturbing symptoms that are unclear to you related to the use of the medical device, a healthcare professional.

This symbol indicates the maximum weight of the user

Manufacturer:

KARMA MEDICAL PRODUCTS CO.,LTD.

Number 2363 section. 2 Da-Shiue Rd., Min-Hsiung Shiang,
Chia-Yi, Taiwan

TEL: 886-(0)5-2066688 FAX: 886-(0)5-2067788

CE Representative

Authorized

representative: 28

Newfield Road, Stoubridge, West Midlands DY9 OJR U.K

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