User manual

MANUEL D'UTILISATION
GEBRUIKSAANWIJZING
GEBRAUCHSANWEISUNG
MANUALE DI ISTRUZIONI
MANUAL DE INSTRUCCIONES
INSTRUKCJA OBSŁUGI
NÁVOD K OBSLUZE

LiNX Memory Positions

For specialised care providers

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Preface

Important note

This manual reflects the latest product developments. Vermeiren has the right to implement changes to this type of product without any obligation to adapt or replace similar products previously delivered.

Pictures of the product are used to clarify the instructions. Details of the depicted product may deviate from your product.

Information available

On our website <http://www.vermeiren.com/> you will always find the most recent version of the following information. Please consult this website regularly for possible updates.

Visually impaired people can download the electronic version of this manual and have it read out by means of a text-to-speech software application.

	User manual of the wheelchair For user and specialist dealer
	User manual of LiNX Memory Positions system For user and specialist dealer
	User manuals of applied operating console and battery charger For user and specialist dealer
	Installation instructions For specialist dealer
	Service manual for wheelchairs For specialist dealer
	EC declaration of conformity

1. Description

The LiNX Memory Positions system allows users and specialised care providers to program preset sitting positions for their wheelchair. These positions can have therapeutic purposes, or they can simply be for the user's comfort.

LiNX Memory Positions uses sensors on several locations on the wheelchair to monitor the sitting position. These sensors ensure that stored sitting positions can be achieved with just one joystick movement. Altering the sitting position one function at a time (backrest, seat, leg rest,...) remains possible.

2. Before use

WARNING

Risk of injuries and damage

- Do NOT install, maintain or use this product unless you have read and fully understood all instructions in this manual and all relevant manuals regarding the wheelchair equipped with this product.
 - Do not modify your product in any way.
 - In case of doubts or questions, do not hesitate to contact your local specialist dealer, your care provider or technical adviser to help you with this.
-  Vermeiren is responsible for software updates/repairs. For adaptations in the software contact Vermeiren, see our website <http://www.vermeiren.com> .

3. Using LiNX Memory Positions

CAUTION

Risk of injuries or damage

- First read the manual of the operating console.

Your wheelchair with LiNX Memory Positions system is delivered with standard configurations. Users and specialised care providers can each program two preset sitting positions to overwrite these standard configurations.

- Comfort (D, C1) : User positions
- Clinical (D, M1) : Therapeutic sitting positions, set by specialised care providers

Therapeutic sitting positions can be used by the wheelchair user, but they can only be altered by the care provider using a dongle to access the settings menu.

-  If the user does not need any therapeutic sitting positions, the specialised care provider can leave these positions unlocked, allowing the user to manage all 4 available positions.


3.1. For users

3.1.1. Using preset positions

CAUTION


Risk of injuries

- All parts of the wheelchair (backrest, seat, leg rest) will move simultaneously once a stored position is chosen. Make sure that no body parts can be trapped or squeezed between moving parts of the wheelchair.
1. Navigate through the S-menu (for adjusting the body position) on the touchscreen until you reach the "Comfort" or "Clinical" submenu.
 2. Press the joystick to the front or back, according to which position you want to assume. The relevant wheelchair parts will start moving.
 3. If you want to stop the automatic movement, shortly press the joystick again. The wheelchair parts will stop moving immediately.

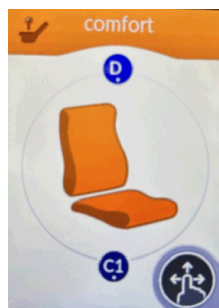3.1.2. Programming preset positions

1. Using the functions in the S-menu, alter the position of the backrest, seat, and leg rest separately until you reach the desired sitting position.
2. When you have reached the body position you want to store in the system, enter the pop-up "Settings" menu (K).
3. Navigate to the "Update position" item. Tap it to open the overview of preset positions.
4. Tap the name of the position you want to overwrite with your new position. Tap "Yes" on the pop-up to confirm the overwrite.
5. You can now use the new preset position anytime by following the instructions in §3.1.1..

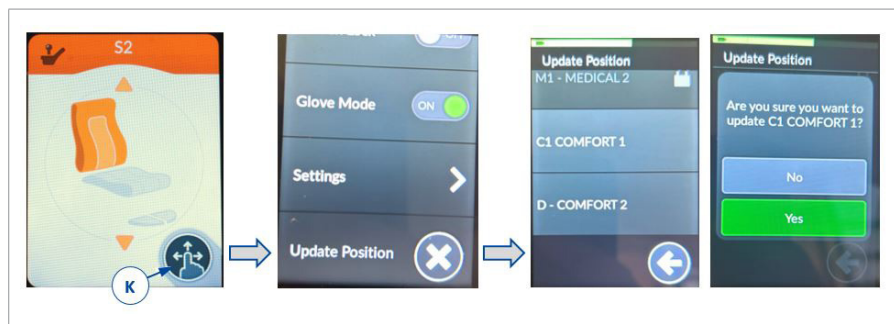
3.2. For specialised care providers

Specialised care providers can preset two therapeutic sitting positions that can be used but not altered by the user. To alter the Clinical positions, a dongle is needed to access the settings menu. Alterations can be done via desktop, or via smartphone.

1. Open the program menu on your desktop or smartphone and go to the "Positions" tab.
2. Open the position you want to update.
3. Under the "Configuration" section of the selected position, set "Allow Occupant to Update Position" from OFF to ON.
4. Upload the program to the wheelchair.
5. Change and save the sitting position of the wheelchair to the desired therapeutic position according to the instruction in §3.1.2..
6. On your desktop or smartphone, set "Allow Occupant to Update Position" back from ON to OFF for the selected position.
7. Upload the program to the wheelchair.

i If locked Clinical positions are not necessary, you can open up all positions to be programmed by the user by following steps 1 -4.

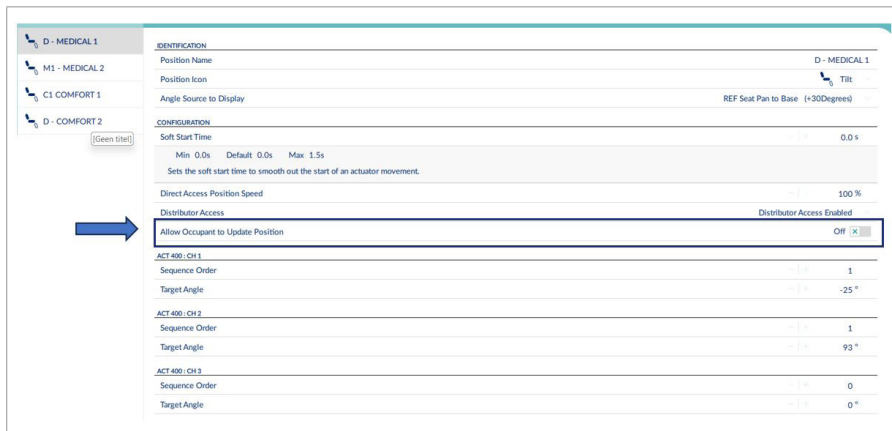
4. Troubleshooting

In case of a possible problem, a system code is displayed by the blinking LED in the ON-OFF button or by a triangle on the touch screen, see the user manual of applied operating console for more information.

If there is an issue with one of the LiNX Memory Positions sensors or its connections, the whole system will be unavailable for use. A large exclamation mark will cover the images in the LiNX Memory Positions menus, and the wheelchair will not react to joystick manipulation.

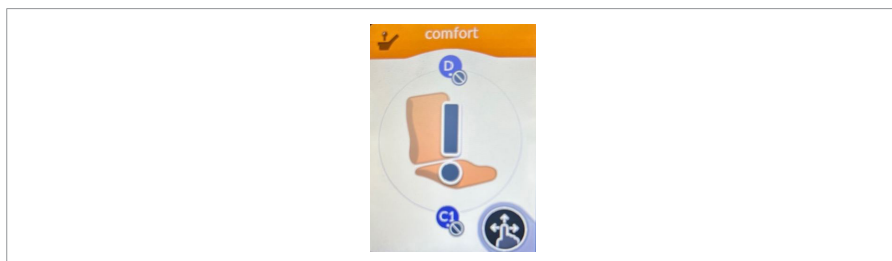Other wheelchair functions remain available as usual.

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SIGMA.SE19.SE39.SE42.SE48

Since 18/01/2024 9:42 (DLX-HKEY01-A I21255047)
1 hours driving0 hours driving / day

Since 18/01/2024 9:42 (DLX-HKEY01-A I21255047)
0 total battery cycles0.1 hours charging / day

No active errors

EVENTS

System Available	Today
Programming Changed	Today
System Available	Today

SYSTEM NAME

SIGMA.SE19.SE39.SE42.SE48

PROGRAM NAME

SIGMA.SE19.SE39.SE42.SE48

CONFIGURATION SUMMARY

6 Drive Functions

12 Seating Functions

1 Mouse Mover Functions

1 Switch Control Functions

6 Motion(s)

4 Position(s)

7 Angle Source(s)

6 Trigger Angle(s)

Gyro supported

Lighting supported

System Summary

D - MEDICAL 1

M1 - MEDICAL 2

C1 COMFORT 1

D - COMFORT 2
(Green title)

IDENTIFICATION

Position NameD - MEDICAL 1

Position Icon

Angle Source to DisplayREF Seat Pan to Base (+30Degrees)

CONFIGURATION

Soft Start Time0.0 s
Min: 0.0s Default: 0.0s Max: 1.5s
Sets the soft start time to smooth out the start of an actuator movement.

Direct Access Position Speed100 %

Distributor AccessDistributor Access Enabled

Allow Occupant to Update PositionOFF

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OEM MVG 2
SIGMA.SE19.SE39....
Flexible Navigation ⓘ v7.0

- Home
- Functions
- User Preferences
- Chair Log
- Modules
- Motions
- Positions**
- Angle Sources
- Trigger Angles
- Core Features
- Drive Limits
- Gyro Limits
- Lighting

SIGMA.SE19.SE39.SE42.SE48

Since 18/01/2024 9:42 (DLX-HKEY01-A I21255047)

1
hours driving

0
hours driving / day

Since 18/01/2024 9:42 (DLX-HKEY01-A I21255047)

0
total battery cycles

0.1
hours charging / day

No active errors

EVENTS	
System Available	Today
Programming Changed	Today
System Available	Today

SYSTEM NAME

SIGMA.SE19.SE39.SE42.SE48

PROGRAM NAME

SIGMA.SE19.SE39.SE42.SE48

CONFIGURATION SUMMARY

- 6 Drive Functions
- 12 Seating Functions
- 1 Mouse Mover Functions
- 1 Switch Control Functions
- 6 Motion(s)
- 4 Position(s)
- 7 Angle Source(s)
- 6 Trigger Angle(s)
- Gyro supported
- Lighting supported

System Summary

- D - MEDICAL 1**
- M1 - MEDICAL 2
- C1 COMFORT 1
- D - COMFORT 2

(Green title)

IDENTIFICATION

Position Name: D - MEDICAL 1

Position Icon:

Angle Source to Display: REF Seat Pan to Base (+30Degrees)

CONFIGURATION

Soft Start Time: 0.0 s

Min: 0.0s Default: 0.0s Max: 1.5s

Sets the soft start time to smooth out the start of an actuator movement.

Direct Access Position Speed: 100 %

Distributor Access: ☐ Distributor Access Enabled

Allow Occupant to Update Position ☐ Off

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All rights reserved, including translation.