User manual

MANUEL D'UTILISATION
GEBRUIKSAANWIJZING
GEBRAUCHSANWEISUNG
MANUALE DI ISTRUZIONI
MANUAL DE INSTRUCCIONES
INSTRUKCJA OBSŁUGI
NÁVOD K OBSLUZE

LiNX Memory Positions

Contents

Preface	4
1. Description	5
2. Before use	5
3. Using LiNX Memory Positions	6
3.1. For users	6
4. Troubleshooting	8

Preface

Important note

This manual reflects the latest product developments. Vermeiren has the right to implement changes to this type of product without any obligation to adapt or replace similar products previously delivered.

Pictures of the product are used to clarify the instructions. Details of the depicted product may deviate from your product.

Information available

On our website <http://www.vermeiren.com/> you will always find the most recent version of the following information. Please consult this website regularly for possible updates.

Visually impaired people can download the electronic version of this manual and have it read out by means of a text-to-speech software application.

	User manual of the wheelchair For user and specialist dealer
	User manual of LiNX Memory Positions system For user and specialist dealer
	User manuals of applied operating console and battery charger For user and specialist dealer
	Installation instructions For specialist dealer
	Service manual for wheelchairs For specialist dealer
	EC declaration of conformity

1. Description

The LiNX Memory Positions system allows users and specialised care providers to program preset sitting positions for their wheelchair. These positions can have therapeutic purposes, or they can simply be for the user's comfort.

LiNX Memory Positions uses sensors on several locations on the wheelchair to monitor the sitting position. These sensors ensure that stored sitting positions can be achieved with just one joystick movement. Altering the sitting position one function at a time (backrest, seat, leg rest,...) remains possible.

2. Before use

WARNING

Risk of injuries and damage

- Do NOT install, maintain or use this product unless you have read and fully understood all instructions in this manual and all relevant manuals regarding the wheelchair equipped with this product.
 - Do not modify your product in any way.
 - In case of doubts or questions, do not hesitate to contact your local specialist dealer, your care provider or technical adviser to help you with this.
-  Vermeiren is responsible for software updates/repairs. For adaptations in the software contact Vermeiren, see our website <http://www.vermeiren.com> .

3. Using LiNX Memory Positions

CAUTION

Risk of injuries or damage

- First read the manual of the operating console.

Your wheelchair with LiNX Memory Positions system is delivered with standard configurations. Users and specialised care providers can each program two preset sitting positions to overwrite these standard configurations.

- Comfort (D, C1) : User positions
- Clinical (D, M1) : Therapeutic sitting positions, set by specialised care providers

Therapeutic sitting positions can be used by the wheelchair user, but they can only be altered by the care provider using a dongle to access the settings menu.


3.1. For users

3.1.1. Using preset positions

CAUTION


Risk of injuries

- All parts of the wheelchair (backrest, seat, leg rest) will move simultaneously once a stored position is chosen. Make sure that no body parts can be trapped or squeezed between moving parts of the wheelchair.
1. Navigate through the S-menu (for adjusting the body position) on the touchscreen until you reach the "Comfort" or "Clinical" submenu.
 2. Press the joystick to the front or back, according to which position you want to assume. The relevant wheelchair parts will start moving.
 3. If you want to stop the automatic movement, shortly press the joystick again. The wheelchair parts will stop moving immediately.

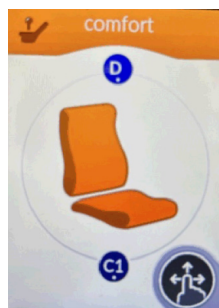3.1.2. Programming preset positions


1. Using the functions in the S-menu, alter the position of the backrest, seat, and leg rest separately until you reach the desired sitting position.
2. When you have reached the body position you want to store in the system, enter the pop-up "Settings" menu **(K)**.
3. Navigate to the "Update position" item. Tap it to open the overview of preset positions.
4. Tap the name of the position you want to overwrite with your new position. Tap "Yes" on the pop-up to confirm the overwrite.
5. You can now use the new preset position anytime by following the instructions in **§3.1.1.**

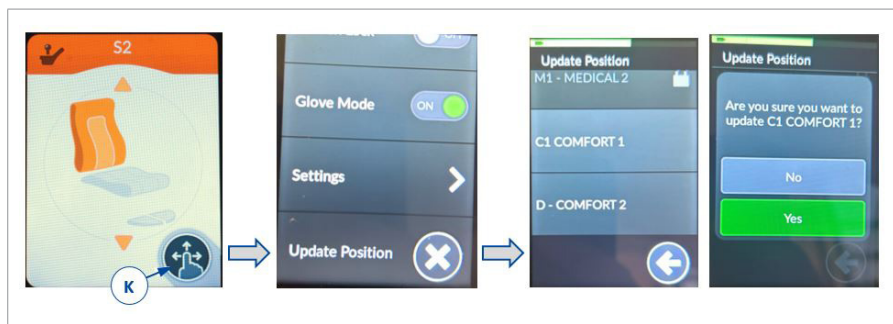4. Troubleshooting

In case of a possible problem, a system code is displayed by the blinking LED in the ON-OFF button or by a triangle on the touch screen, see the user manual of applied operating console for more information.

If there is an issue with one of the LiNX Memory Positions sensors or its connections, the whole system will be unavailable for use. A large exclamation mark will cover the images in the LiNX Memory Positions menus, and the wheelchair will not react to joystick manipulation.

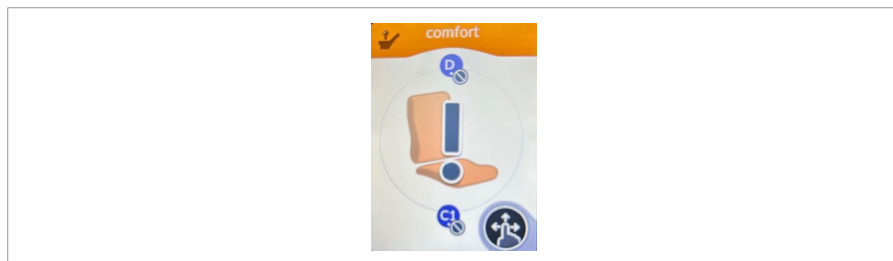Other wheelchair functions remain available as usual.

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All rights reserved, including translation.