Installation manual

MANUEL D'INSTALLATION
INSTALLATIEHANDLEIDING
INSTALLATIONSANLEITUNG
MANUALE DI INSTALLAZIONE
MANUAL DE INSTALACIÓN
INSTRUKCJA INSTALACJI
PŘÍRUČKA PRO INSTALACI

Blumil GO, Blumil Sport

Preface	4
1. Scope of delivery	5
2. Installation	5
2.1. Installing the mounting bar receivers	5
2.2. Adjusting the Blumil receiver to your wheelchair	8
2.3. Brake adapter plate	11

Preface

Important note

Pictures of the product are used to clarify the instructions. Details of the depicted product may deviate from your product.

Information available

On our website <http://www.vermeiren.com/> you will always find the most recent version of the information in this manual. Please consult this website regularly for possible updates.

Visually impaired people can download the electronic version of this manual and have it read out by means of a text-to-speech software application.

	User manual For user and specialist dealer
	User manual of the wheelchair For user and specialist dealer
	Installation instructions For specialist dealer
	Drawings of (spare) parts For specialist dealer
	EC declaration of conformity

1. Scope of delivery

The following components are part of the delivery:

- Blumil device for foldable wheelchairs
- Mounting bar
- Manual

Before use, check if everything is included and that no products are damaged.

Please note that the basic configuration may differ in different European countries. Contact your specialist dealer.

2. Installation

Blumil can be connected to a foldable wheelchair using the removable mounting bar. The only parts that are permanently added to the wheelchair are the receivers for the mounting bar.


2.1. Installing the mounting bar receivers

- i** Different types of receivers are available, depending on the wheelchair frame. Make sure to use the correct type for your wheelchair.


2.1.1. Trigo wheelchairs

i Applicable to:

- Trigo S
- Trigo T

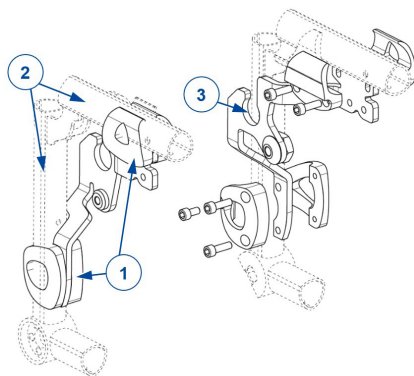
1. Place the clamps (1) of the receiver over the frame tubes (2) of the wheelchair (one below the seat, one to the front frame tube).
Ensure that the receiver is bending inwards. The slit (3) for the mounting bar needs to be oriented upwards.
2. Slightly tighten all screws on the clamps so the receiver remains in place.
3. Repeat for the other side.

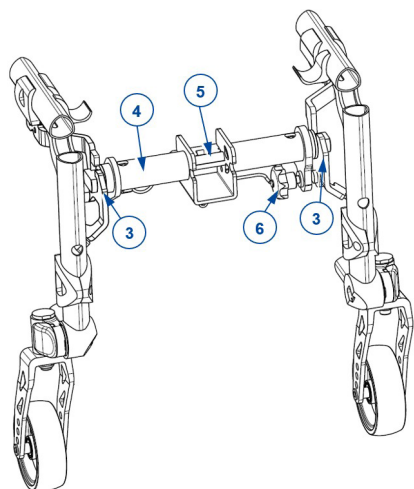
4. Insert the mounting bar (4) into the slits (3), and rotate to click into place. Make sure that the bar is oriented with the Blumil connection point (5) to the front.
5. Tighten the star knob (6) by hand.
6. Firmly tighten all screws to secure the clamps into place.

2.1.2. V-series

- i** Applicable to:
- V500 Light
 - V500 Light 30°
 - V500 Active

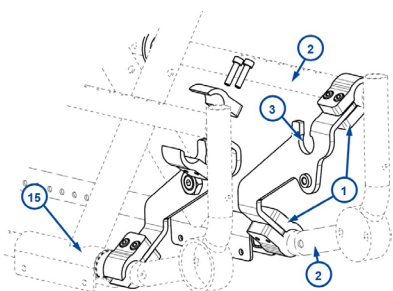1. Place the clamps (1) of the receiver over the frame tubes (2) of the wheelchair (one below the seat, one to the lower frame tube). Ensure that the lower clamp is placed against the plastic cross attachment (15) on the lower frame tube.
Ensure that the receiver is bending inwards. The slit (3) for the mounting bar needs to be oriented upwards.
2. Slightly tighten all screws on the clamps so the receiver remains in place.
3. Repeat for the other side.

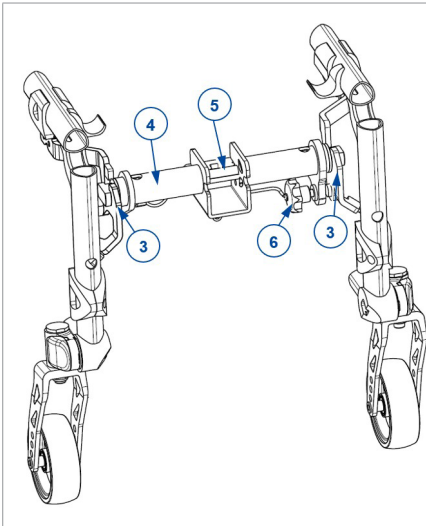
4. Insert the mounting bar (4) into the slits (3), and rotate to click into place. Make sure that the bar is oriented with the Blumil connection point (5) to the front.
5. Tighten the star knob (6) by hand.
6. Firmly tighten all screws to secure the clamps into place.


2.2. Adjusting the Blumil receiver to your wheelchair

- i** Blumil Sport: Make sure to mount the standing support before adjustment.


To adjust the height and angle of the Blumil connection:

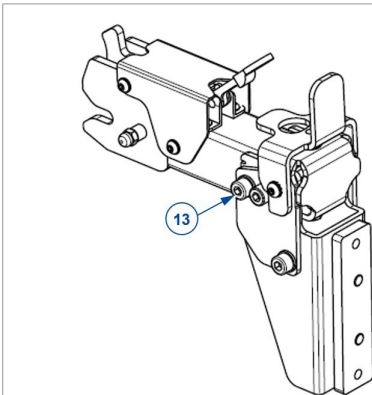
1. Slightly loosen the bolt-nut connection (13) to adjust the Blumil connection.

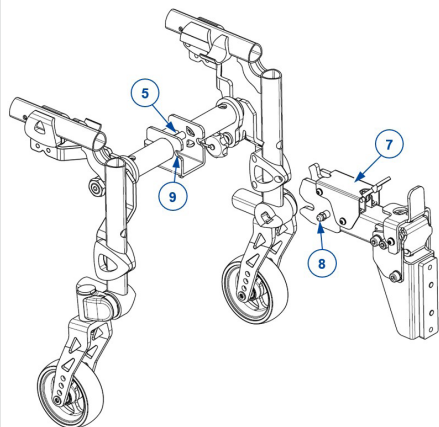
2. Slide the Blumil connection (7) to a height that corresponds to the Blumil connection point (5) on the mounting bar.
3. Try to click the connection into the connection point to check the angle. You should be able to smoothly click the pin (8) of the Blumil connection into the slot (9) of the Blumil connection point.

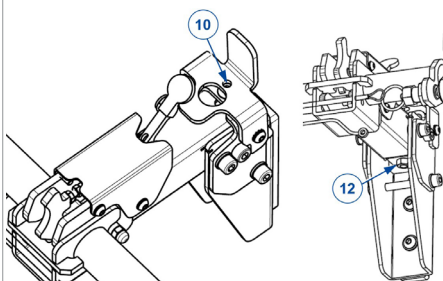
4. If necessary, adjust the angle of the connection by turning the top screw (10). Make sure that the nut (12) on the bottom of the screw is loosened first.

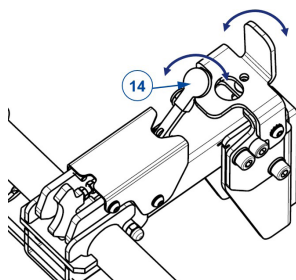
5. Once the connection smoothly clicks into the connection point, try out the locking lever (14) to ensure that it works well.
6. If the red safety lever easily opens again with a flick of your thumb, and the connection easily slides back out of the connection point, you have determined the correct settings.

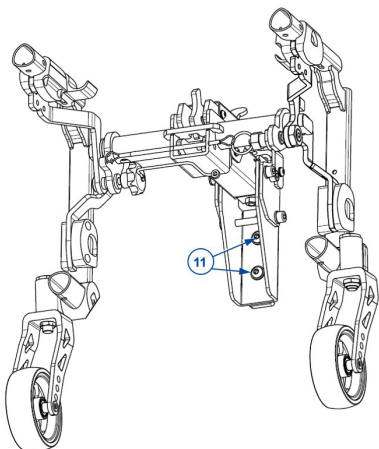7. With Blumil still connected to the wheelchair, firmly secure the height adjustment by tightening the 2 screws (11) on the inner side of the Blumil connection.
8. Remove Blumil a bit from the wheelchair again.

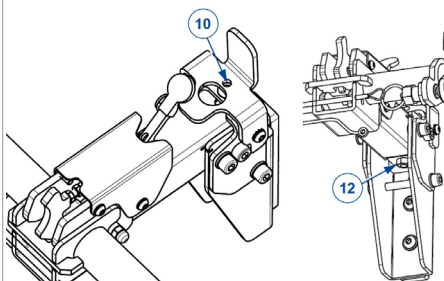9. Firmly secure the angle adjustment by tightening the nut (12) on the bottom of screw (10). Use the Allen key to prevent the screw from turning while tightening the nut.

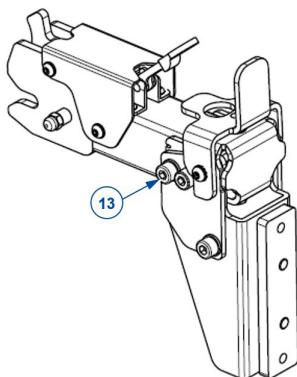10. Firmly tighten the bolt-nut connection (13) to secure the whole adjustment system.


2.3. Brake adapter plate

i Applicable to:

- Trigo S
- Trigo T

If the Blumil receiver interferes with the standard brake position, the brakes need to be repositioned to the optional brake adapter plates.

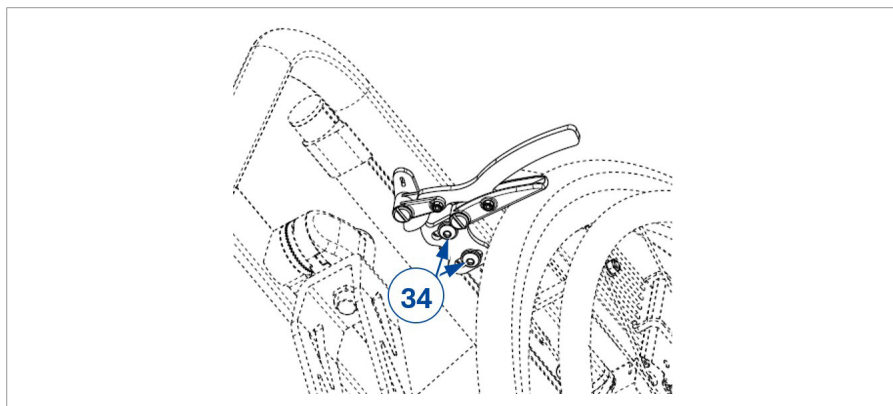- Adapter plate A can be taken as standard.
 - Adapter plate B can be used in configurations where adapter plate A is insufficient for correct brake placement due to rear wheel camber, rear wheel position, ...
1. Remove both screws (34) and take the brake with attachment off the wheelchair frame.
 2. Mount the receiver (35) as described in §2.1.1..
 3. Mount the brake adapter plate (36) to the receiver by inserting the bolt-nut connections (37). Tighten firmly by hand.
 4. Mount the brake into the slot (38) of the brake adapter plate. Tighten the screws (34) firmly by hand.
 5. If necessary, adjust the brake by sliding it forward/backward in the slot.
 6. Repeat for the other side.

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All rights reserved, including translation.